
SUPER NOVA

HIGH FREQUENCY MASSAGER

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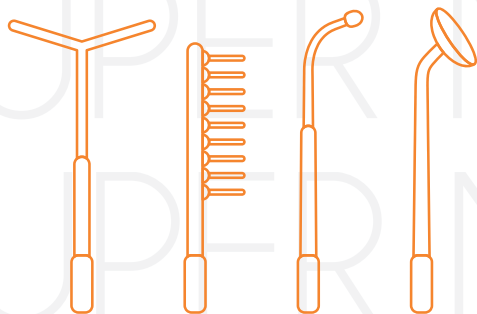USER'S MANUAL

Portable professional
high frequency device

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GESS[®]
Expert

We thank You for purchasing the massager «Super Nova»!

This device is not a medical equipment and is to be used for prophylactic purposes.

We reserve the right to modify the design and make amends to the description of our products. Changes can be carried out without prior notice, in regards to its color and components.

Please, read the instruction manual before use. After reading, keep the manual in a good state for further reference to it if required.

Pay particular attention to the section marked «**Safety precautions**».

English

Deutsch

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Content

Safety precautions.....	3
What is Darsonvalization and how does it work?.....	4
Description of the device	6
Conducting Darsonvalisation procedures.....	7
Possible malfunctions and their troubleshooting	10
Technical specifications	10

Safety precautions

- Carefully read the instruction manual before using the device!
- Do not run the device in rooms with high humidity. It is not recommended to use the device in bathrooms and swimming pools.
- After moving the device to a warm room from a cold one, it is only allowed to be used not earlier than 1 hour.
- Before connecting the device to the mains, make sure that the glass electrodes, the encasing of the device and the power cord insulation are not damaged. If any damage is found further operation of the device is prohibited.
- It is prohibited to turn the device on if no electrode is connected to it. Do not operate the device if its casing gets into contact with water or other liquids. Any liquid from the device casing should be carefully removed with a soft cloth or tissue.
- It is strictly prohibited to use the device if water (or any other liquid) gets inside the device casing. In this case, You should contact the service center for inspection and repair.
- It is prohibited to touch the surface of the electrodes in the area containing the metal holders on the device-casing when the device is operating.
- Any treatment for children should be carried out by a competent adult person after consulting a doctor.
- The device should be stored out of the reach of children. Do not leave the device unattended!
- It is recommended to keep the device away from direct sunlight and away from places with high temperature or humidity.
- Protect the device, and especially the glass electrodes from impact.

What is high frequency massager and how does it work?

Darsonvalization is based on the therapeutic effect of an alternating current with a high voltage on the human body, through a gas-filled glass electrode.

The effect is caused by pulsed alternating currents at high frequency (100-100 kHz) and high voltage (about 2 kW), but at a relatively low amperage (up to 100-200 mA). Alternating currents at low amperage do not lead to ionic changes in the membranes, do not lead to increases in temperatures and do not cause muscle tissue to contract. The positive effect is achieved by improving the outflow of venous blood and lymph, thus accelerating metabolism in tissues.

Using the Darsonval pulse massager activates biochemical metabolism in the skin and subcutaneous layers, normalizes blood circulation, improves tissue tropism and supply them with oxygen, reduces the threshold of sensitivity of nociceptors to external stimuli, providing an analgesic effect.

Due to the activation of metabolic processes and blood circulation, the devices improves the penetration of therapeutic and cosmetic products into the skin and the subcutaneous layers and improves muscle tone.

Generalized darsonvalization helps to tone the brain vessels, eliminates headaches, normalizes sleep and increases work efficiency.

There are two ways to use the Darsonval impulse massager at home:

What is high frequency massager and how does it work?

- Contact («silent discharge») - the vacuum electrode of the massager has an effect on a particular area of the skin, by being in direct contact with it. Treatments in the «silent discharge» mode normalize sebum secretion, improve skin elasticity and turgor, prevent the development of wrinkles.
- Contactless («spark discharge») - the massage nozzle is located at a distance of 1-4 mm from the skin surface. When the electrode is taken off the skin, a spark discharge is produced, the air around the massager gets filled with ozone, and a slight tingling sensation is felt on the treated area of the skin.

The first method is applied for the treatment of aging and problematic skin, cellulite, varicose veins, and to improve hair growth, and the second - for pain relief, treatment of diseases with a neuralgic character, resorption of hematomas and promote wound healing. In addition, both the spark and silent darsonval discharge have a pronounced bactericidal effect. Desired effect is achieved in both cases, on average after 10-12 treatments.

Regular use of this device improves vascular tone and strengthens the immune system, normalizes the central nervous system: improves sleep, increases efficiency, reduces headaches and fatigue.

INDICATIONS FOR DARSONVALIZATION.

This electrotherapy method is widely used in medicine, that is why, the use of Darsonval pulse massager has quite an impressive list of indications

A partial list of indications for the application of the Darsonval device :

1. **Cosmetology and dermatology**

Effects on the skin of the face, neck and body:

- pustules, impure skin, inflammatory infiltrates acne
- fat porous skin
- mature, flabby, wrinkled skin with reduced turgor
- pale lifeless skin with reduced nutrition
- psoriasis
- Lichen planus
- eczema
- erythema
- pruritus, including allergic ones
- red grainy nose
- keloids
- focal and systemic scleroderma
- herpetic eruptions
- keloid and postoperative scars
- postoperative edema and infiltration
- neurodermatitis

Effects on the scalp:

- intensive hair loss, seborrhea, dandruff
- seborrheic alopecia and other types of alopecia
- strengthening the hair follicles

What is high frequency massager and how does it work?

- stimulation of hair growth
- Effects on areas of the body:
- all types and stages of cellulite.

2. Combined therapy

- Diseases of the peripheral nervous system, including those with pain syndrome, neuroses: sciatica, neuromyositis, neuralgia, osteochondrosis, climacteric neurosis, migraine, sleep disorders, myopathy, inflammation of the nerve plexus, neurasthenia, Raynaud's disease, acroparesthesia, cerebral atherosclerosis, effects of spinal cord and brain trauma, antipruritic effect for dermatitis and allergies.
- Peripheral vascular diseases: sores and other trophic disorders, varicose veins of the lower extremities.
- Joint diseases: polyarthritis arthritis, spondylarthritis.
- Diseases of the nose: sinusitis, allergic rhinitis.
- Diseases of the respiratory system: bronchitis, tracheitis, asthma.
- Diseases of tendons, joint capsules, muscles and aponeuroses: myalgia, myositis.
- Gynecology: colpitis, vaginismus, erosion of the uterus, frigidity, itching of the vulva, endocervicitis, postpartum perineal rupture, cracked nipples of mammary glands.
- Urology: impotence, urinary incontinence, cystitis, prostatitis.
- Dentistry: pulpitis, glossalgia, gingivitis, parodontopathy, stomatitis.
- Non-healing wounds.

CONTRAINDICATIONS FOR USING HIGH FREQUENCY CURRENTS.


The use of darsonvalization is not recommended in the following cases:

- malignant neoplasms
- benign neoplasms with a tendency to progression
- bleeding
- presence of a pacemaker
- blood clotting disorders
- skin sensitivity disorders
- epilepsy
- thrombophlebitis
- pregnancy
- cardiac arrhythmia
- expressed forms of rosacea (for example a red dilated net of vessels on the face)
- hirsutism - excessive hair growth on the face
- active pulmonary tuberculosis
- fever
- cardiovascular insufficiency
- purulent processes before surgery
- individual intolerance to electric current
- alcohol intoxication

What is high frequency massager and how does it work?

CAUTION! Despite the fact that each darsonval device comes with a set of instructions, supplied with a detailed list of indications and contraindications and describing all the ways and methods the device affects the body, for medical purposes, it is required to have a consultation with the attending physician about the device prior to darsonvalization.

Description of the device

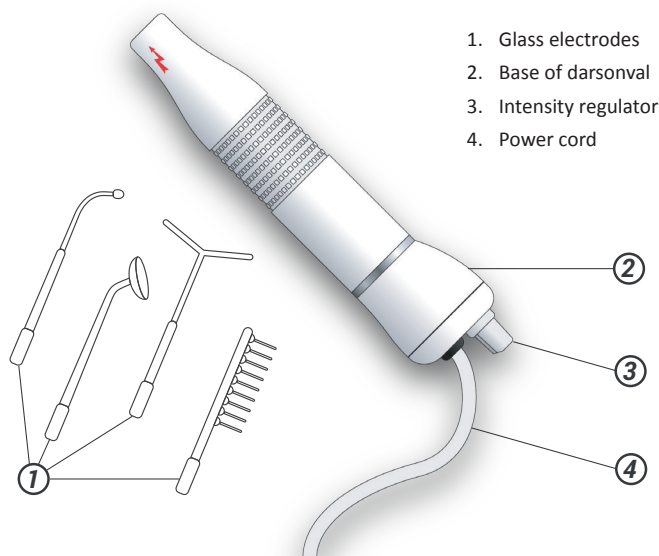"Drop" Electrode – designed for point treatment, used for the treatment of local inflammation of the skin.

"Fork" electrode – ideal for legs, arms and neck treatment, used for the prophylactics of varicose, cellulite and osteochondrosis.

"Mushroom" electrode – designed for generalized procedures for skin, body and face treatment.

"Comb" electrode – used for hair and scalp treatment.

Conducting Darsonvalization procedures

RULES FOR PULSE MASSAGER APPLICATION

1. Select the desired electrode and gently insert it into the metal holder until it fits perfectly. The electrode should be installed and removed with a slight effort. Do not use force! This can cause damage to the electrode.
2. Before connecting the device to the mains, make sure that the control knob is turned to the maximum anti-clockwise until it clicks into the «OFF» position.
3. Switch the device on by turning the knob a little clockwise (until it clicks). Bring the electrode tip gently against the skin surface. Set to the desired effect force by smoothly turning the knob clockwise - the extreme clockwise position of the device knob corresponds to maximal intensity of exposure. Be careful!
4. Caution! It is necessary to switch the device off before changing the electrode by turning the intensity regulator to the maximum anti-clockwise position of the arrow until it clicks.
5. Service the device only after completely disconnecting it from the mains! Dirt from the surface of the device and glass electrodes must be removed with a soft damp cloth. It is allowed to use a neutral detergent or alcohol solution to clean the electrodes. The casing of the device and the electrodes must be thoroughly dried before further use.

METHOD TO CONDUCT

DARSONVALIZATION PROCEDURES

When conducting a darsonvalization procedure, the skin should be clean and free of makeup, creams and others. Therefore, before the darsonvalization procedure it is recommended to thoroughly clean the skin.

Exposure is achieved with slow light movements about the main massage lines of the face and body. Exposure in the «silent discharge» mode is carried out by direct contact - no air gap between the electrode and the skin surface. When exposed in this manner sensations may be weak, and a pleasant warmth may be felt.

Exposure in the «spark discharge» mode is achieved with an air gap of 1-4 mm between the electrode surface and the skin. When exposed in the «spark discharge» mode, sensations can be as follows: tingling, silent crackling, intense ozone smell. The blue sparks can be noticed, resembling thin little lightning bolts.

It is recommended to sprinkle talcum or baby powder on the body before the procedure to improve the sliding of the electrode on the skin, for direct exposure. When exposed to the face - either cosmetic talcum or powder. When exposed to wounds and ulcers - through a clean, dry gauze. In order to improve the supply of the skin with oxygen special procedures can be used by applying an oxygen rich gel onto the skin. Exposure is started at low intensity and is gradually increased to the desired value.

Once completed, turn off the device by turning the knob anti-clockwise until it clicks and disconnect it from the mains. Service the device and store it in a storage case.

CONDUCTING COSMETIC PROCEDURES

Procedures on the face and neck

1. Remove all the makeup, eliminating any cosmetic product from the skin using a milk or tonic for Your skin type. If suffering from acne and oily skin, a deeper cleansing is recommended using an exfoliating scrub or peeling with a special sponge.
2. Remove excess moisture from the skin tissue and apply a little cosmetic talcum to improve the sliding of the electrode (if suffering from

Conducting Darsonvalization procedures

3. acne, it is better to use talcum with anti-inflammatory drying components).
4. Improving skin nutrition. Exposure is achieved using a «mushroom» glass electrode, and is performed with slow sliding movements, in an upward direction - from the chin to the earlobe, from the nose to the temples, from the eyebrows up to the scalp, i.e. in accordance with the basic massage lines. In the neck area, movements are carried out on the massage lines from bottom to top, without affecting the area of the thyroid gland. The intensity of exposure is as follows: To activate blood circulation while reducing skin turgor, on skin with wrinkles, low nutrition and saggy skin - using the method of «silent discharge» with weak or moderate intensity. To normalize the activity of sebaceous glands and reduce the porosity of the skin - using the method of «silent discharge» with moderate or high intensity.
5. Treatment of problematic skin. When exposed to acne, pustular rash, inflamed scleroderma, a «Drop» electrode should be used. Exposure is directly on the problem areas and slow movements are made on the area around it. The intensity of exposure is as follows: When exposed to purulent infiltration - using the «silent discharge» (contact) method, with an average intensity. When exposed to pustules - using the «spark discharge» method with an air gap of 1-4 mm at a high intensity to have a caustic effect.
6. After the procedure is completed, it is recommended to apply a nourishing moisturizer corresponding to Your skin type to the face and neck, you can also perform a manual massage of the face and apply a nourishing mask for 15-20 minutes. Darsonvalization performing procedures can be combined with ultrasound massage (the same day) and electro-therapy (on different days).
7. Procedures are conducted daily or every other day. The total duration of one session is 5-10 min (the duration should be gradually increased). The duration of exposure to the neck region - should not be more than 5 minutes. One course may consist of 10 procedures - to improve the nutrition of the skin, and 20 procedures - to treat problematic skin. Treatment courses can be repeated 3-4 times a year.

Effect on the scalp

(the stimulation of hair growth and the treatment of migraine)

1. Hair must be clean and dry. It is imperative to remove any hair barrettes and comb the hair well. Exposure is achieved using the «Comb» electrode and the contact method (without any air gap) at an average intensity.
2. The electrode must be smoothly and slowly moved from the front to the back, and if the hair is short, then it should be moved in the opposite direction. The total procedure duration should be 10-15 minutes.
3. It is very effective to rub a therapeutic serum into the scalp or make a hair mask when the exposure is completed.
4. Procedures should be carried out daily. The course of the treatment is 20-30 procedures. The course can be repeated 3-4 times a year.


Treating varicose veins, vasomotor spasm

Exposure is achieved using the «fork» electrode and the contact method at an average intensity. It is recommended to sprinkle talcum powder on the area to be treated.

1. Exposure is achieved with slow movements from the tips of the toes up to the knee for a duration of 5-10 minutes, then from the heel to the popliteal fossa for a duration of 3-5 min., then from the knee to the groin for a duration of 5-10 min., and finally from the popliteal fossa to the gluteal folds for 5-10 min. The movement patterns are shown in the figures.
2. When treating varicose veins, after the procedure, it is recommended to apply any gel or cream used to treat varicose veins.

Conducting Darsonvalization procedures

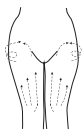- It is recommended to perform the procedures daily or every other day. The course consists of 15 procedures.

Treatment of cellulite

Exposure is achieved using the «mushroom» electrode and the contact method on the area with cellulite and dilation of superficial dermal capillaries. The exposure intensity is selected from medium to high. After the procedure it is recommended to apply anti-cellulite gel or cream.


Exposure to the hands

The procedure is performed using the «mushroom» electrode and the contact method (without air gap) at an average intensity. Exposure is achieved through smooth circular motion in a spiral from the fingers to the shoulder area. The procedure is particularly effective for the treatment of Raynaud's syndrome. Typical schemes are given as shown in the figure.

Exposure to trophic ulcers and non-healing wounds

It is recommended to perform the procedure when changing the wound dressing.

- Cover the surface of the wound or ulcer with a clean dry gauze. Exposure is achieved using the «mushroom» electrode and the «spark discharge» mode with an air gap of 1-4 mm., at an average intensity. Movements are carried out along the edge of the wound or ulcer with a radius of 5-7 cm. Exposure duration amounts to 10 min. The course consists of 10-15 procedures.

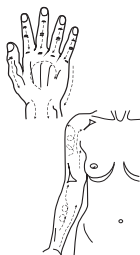- After the procedure on the wound or trophic ulcer area, apply the medication recommended by your doctor (ointments or gels).

Exposure to the collar region

(when treating neuroses, sleep disorders)

The procedures are performed under the doctor's directions. Exposure is achieved using the «mushroom» electrode at the «silent discharge» mode and at low intensity. Movements are performed on the neck region from the back of the neck to the shoulders. Exposure duration amounts to 5-10 minutes. Procedures are recommended daily or every other day. The course consists of 10-15 procedures.

Effect on pain area

(myositis, neuritis, neuralgia, sciatica).

Procedure is recommended after a doctor's consultation. Exposure is carried out directly on the painful area at the «spark discharge» mode with a gap of 1-4 mm. The intensity of the exposure is average. Treatment duration amounts to 10 min. The course consists of 10-15 procedures.

NOTE: In all cases of darsonvalization procedures for medical purposes, you must first consult with your doctor!

Possible malfunctions and troubleshooting

Malfunction	Possible causes of the malfunction	Troubleshooting
The device is not functioning	The device is not connected to the mains	Connect the device to the mains, check the reliability of the plug in the socket.
	The intensity regulator is set to «OFF».	Turn the knob of the regulator clockwise until it clicks.
	The leakproofness of the gas-filled glass electrode was compromised	Try setting another electrode. If necessary, contact your service center or purchase a new electrode. Warning! Never touch the jack when setting the electrodes when the device is ON.
	The device is faulty	Contact the service center for inspection and repair.
The effect of the device became weaker, discharge more «scratchy»	The surface of the glass electrode is dirty	Clean the electrode as recommended in the «Working rules»

Technical specifications

Model	GESS-623
Signal form	pulse with a high-frequency filling
Pulse repetition rate	100 Hz $\pm$ 10 Hz
Frequency of the RF filling	110 $\pm$ 50 kHz
Average output current, does not exceed	150 mA
Output voltage	not more than 8 kW, up to max of 16 to 25 kW
Power source	mains 220 V $\pm$ 22 V, 50 Hz $\pm$ 0,5 Hz
Power consumption	60 W
Packing size	280x60x232 mm
Packing weight	0.78 kg

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Room 101, Pudong New Area District, Huanlindun street, side street 799, № 12, Shanghai, China.

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Authorized organization by the manufacturer in Europe:

GESS MARKET OÜ

Reg.code: 12706981; VAT No.: EE101758947

Room No. 219, 2nd floor, Tulika 19 (A) 10613 Tallinn, Estonia | EU

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WhatsApp: +7 (911) 160-11-42

Warranty Maintenance
(provided warranty provisions
are complied with) includes
warranty repair of the items
within the warranty period
of 1 year or replacing them
with equivalent good ones if
repairing is impossible.

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