Rowenta

FR

EN

ES

DE

IT

PT

UK

CS

NL

PL

RU


SK

RO

HU

YOUTH FACIAL MASSAGER

DESCRIPTION

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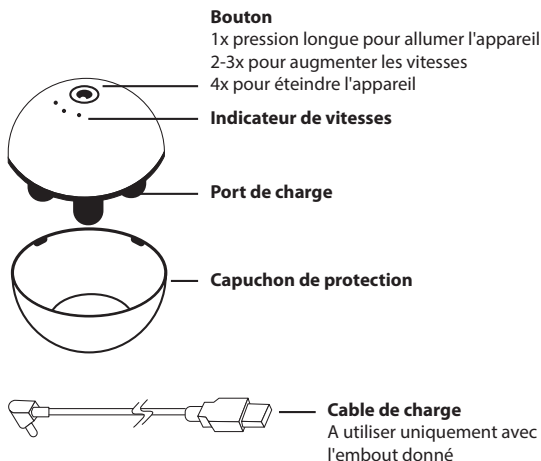
DESCRIPTION

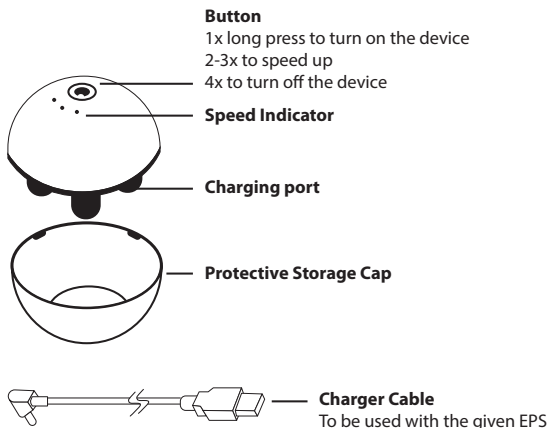BEFORE FIRST USE

- You have purchased the Rowenta Youth Facial Massager: Thanks to its microstimulation, this appliance will restore the glow of your skin and slow down the first signs of ageing. We advise to use it 2min/day 4times/week, in the morning or at night, after applying your usual cream to get the best results.
- To start, press the button. Personalize your routine choosing the speed that you prefer : 1st speed is already set off when you start the device, press 1-2 more times to speed up. Press the button one more time to turn off the device. Use the Youth Facial Massager avoiding sensitive areas (lips, eyes) and creating mouvement from the bottom up.

SAFETY INSTRUCTIONS

- For your safety, this appliance complies with the applicable standards and regulations (Low Voltage Directive, Electromagnetic Compatibility, Environmental...).
- As an additional safety measure, we recommend the installation of a residual current difference device in the electrical circuit supplying the bathroom, acting at not more than 30 mA. Ask your installer for advice.
- For countries subject to EU regulations (CE marking): This device can't be used by children. This appliance can be used by persons with impaired physical, sensory or mental abilities, or lacking experience or knowledge, if they are properly supervised or if they have received instructions concerning how to safely use the appliance and if they have understood the risks involved. Children should not be allowed to play with the appliance. Children should not clean or perform maintenance procedures on the appliance. The battery must be removed from the appliance before it is scrapped.

The appliance must be disconnected from the supply mains when removing

the battery. The battery is to be disposed of safely.

- For other countries not subject to EU regulations: This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance. This appliance contains batteries that are only replaceable by skilled people.
- Your device should be used only for domestic use. It could not be used for professional purpose. Warranty become null and void if used wrongly.
- Before any use, please ensure that the product and accessories have any defect. If your device is damaged or shows any defect, don't use it anymore. If your appliance fall down and works abnormally after, don't use it and contact after sales service. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.

CAUTION : Keep your appliance in a dry place. Never use your appliance with wet hands or near water contained in bath tubs, showers, sinks or other containers... Do not

immerse or put under running water the appliance even for cleaning purposes. Never use your device with wet hands.

- You can however wash the fingers under water after use. Make sure to dry it with a towel after all to remove all stain of water.
- Never use the dishwasher to clean the appliance. Don't clean your appliance with abrasive or corrosive products.
- Don't use this appliance when the temperature is below 15°C and above 35°C. Always keep it on a flat surface resistant to heat. Never put it in contact with other inflammable materials. If the appliance becomes very hot during use, avoid skin contact.
- To avoid heating, you should use the appliance slowly and without effort.
- Detachable supply unit

WARNING: For the purpose of recharging the battery, only use the supply unit and base 0952 provided with this appliance. For an optimal charge, the device should be correctly inserted into the base. Make sure your appliance is in accordance with the supply network of your country. Don't unplug the cord by cord pulling but by plug pulling. Don't use Electrical extension cord.

- The appliance must be unplugged after battery full charge, before cleaning or in case of abnormal functioning.
- This appliance is a skincare appliance, not a medical device. If you are on medication or if you suffer from medical problems, please contact your doctor before using our appliance. Make sure to use the device only on face skin, not on eyes or body. Always use the device in motion on your skin.
- Do not use the appliance:
 - With exfoliating cosmetics, peelings, homemade cosmetics, essential oils, medical creams and corrosive products
 - on wet, chapped skin, open wounds, or when you suffer from eye diseases, eye irritation, eye infection, etc...
 - on the body, directly on eyes, eyelashes and eyelids. Its use is exclusively reserved for the face area
 - on damaged and sun burned skin
 - on unhealthy skin
 - on areas which may present a health hazard (varicose veins, moles, birthmarks)
 - if you suffer from skin conditions and diseases such as eczema, in-

flammation, acne...

- If you have done a cosmetic surgery on your face (botox, lifting, etc...)
 - If you have a pacemaker or dental implant
 - If you are allergic to some metals (nickel)
 - If you take steroids or photosensitive substances
 - If you are pregnant
 - If you have any tattoo on your face or permanent makeup
- Please stop using the appliance in case of allergy, irritation, redness or if you feel pain or any signs of discomfort. If you wear piercings on your face, please take it off before using the device. Take care to not use the device next to those areas. Before use make sure to take glasses, contact lenses and jewels off too. We advise to tie your hair before any use..

CLEANING / RECYCLING

ENVIRONMENT PROTECTION FIRST :

Your appliance contains valuable materials which can be recovered or recycled. Leave it at a local civic waste collection point.

Recycling the batteries:

Before removing the LI-ION batteries, discharge them completely. Open the battery compartment using a tool, cut all the connections and then recycle the batteries in a safe manner.

Those instructions are also available on our website
www.rowenta.com

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