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# IMPORTANT SAFEGUARDS

## Statement

First time use may accompany itch or slight bleeding due to the fact that your gum just starts getting used to the high frequency cleaning of electric toothbrush. Kindly suggest begin with Soft Mode or Sensitive Brush Head. Generally speaking, this phenomenon will disappear after 1-2 weeks. If you do not feel any significant improvement after more than 2 weeks, please stop using this product and promptly consult your dentist. This product is not suitable for users that under the age of 6.

## \* Tips

1. Do not share brush head with others.
2. Do not wrap the charging cable around the charging base while charging.
3. Do not put the travel case in water or any liquid.

## \* Notes


1. This product is designed for oral cleansing. Do not use it for any other purpose. If you experience any discomfort or pain, please stop using the product and consult your dentist.
2. Make sure to check this product for any apparent damage before using it. If damaged or dropped, the product should not be used.
3. Before connecting to the power supply, please check whether the local power supply voltage meets the specification of the product.
4. Do not use the product if any damage is found on the charging base.
5. Replace brush head every 3 months.

## \* Warnings

1. Consult your dental professional before use if you once had orthodontic treatment, dental prosthesis or implantation.
2. This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer before use.
3. This product contains a built-in rechargeable battery. Please keep the product away from flames, and avoid charging, using, or placing it in high temperature.
4. Don't put the charging base in water or other liquid.
5. Make sure the power plug and interface are completely dry before connecting to the power supply.
6. Do not bump or drop the product. If it is damaged due a collision or being dropped, please stop using the product.
7. Do not modify, disassemble or repair the product by yourself.
8. Do not use the product in shower or bath.
9. Do not place metal objects on the charging base.


# I. Introduction

Thank you for choosing DR.BEI to assist on your oral health. Please see the list below to check if you have all the accessories.

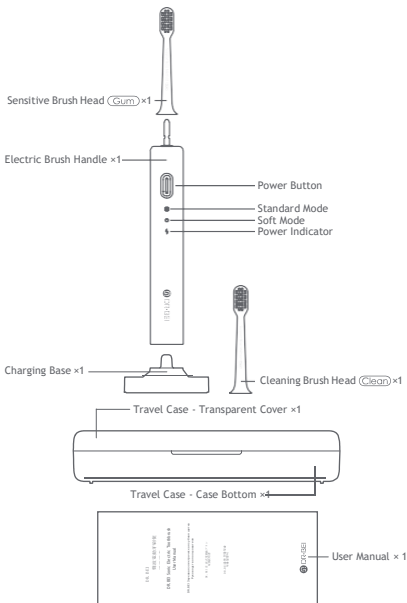
DE.BEI provides you two brush heads for different oral conditions.

**Sensitive Brush Head (Gum)**

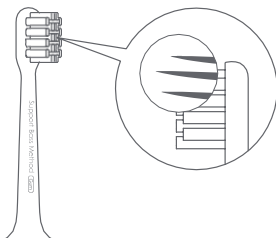**Kind Reminder \***: Sensitive brush head for people with oral sensitivity and gingival bleeding.

Cleaning Brush Head (Clean)

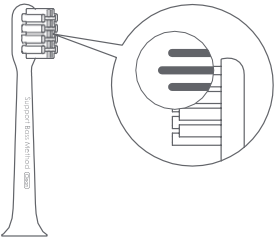Kind Reminder \*: Cleaning brush head is suitable for people with more teeth loss, such as cigarette stains and tea stains.

\* Above instructions are only for your reference. Please choose brush head that brings you the most comfort as oral conditions differs a lot from people to people.

## II. How to Use

### Tip


First time use may accompany itch or slight bleeding due to the fact that your gum just starts getting used to the high frequency cleaning of electric toothbrush. Kindly suggest begin with Soft Mode or Sensitive Brush Head and changes later on your needs.

### a. Install Brush Head

Please follow the picture below:

Ensure the bristle point to the heading direction of the button on the handle and then proceed with installation.

Kindly advise that you pick up the handle from the groove beside the power button and remove it from the groove as the bottom of the handle is fixed in travel case.

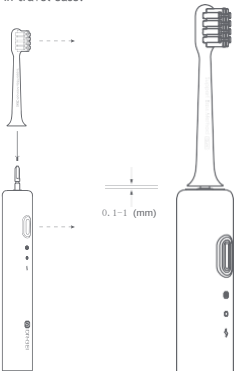
b. Modes Switch

Turn On: Press power button. Default Standard Mode for first boot.

Modes Switch: Press power button again to change modes 5 seconds after power on.

Memory Mode: The system will record the last mode used for more than 5 seconds as Memory Modes and turn to it next time.

Turn Off: Press power button to turn off after at least 5 seconds use or automatically turned off after 2 minutes.

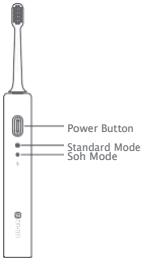c. How to Brush Teeth


DR.BEI's Sonic Electric Toothbrush enables a better support for Bass method recommended by American Dental Association (ADA) and effective cleaning of plaque in your mouth.

The following suggestions can help to deliver superior cleaning results :

Kindly advise that you press the brush against your teeth with a force less than 2N and move at every 2 teeth in your first use. With the 30s smart reminder, you will get a thorough cleaning of both inside and outside of your teeth.

d. Charging

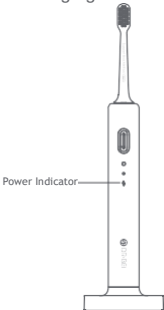Power Indicator Green: 20%-100% power left, feel free to use.

Power Indicator Orange: 10%-20% power left, please charge as soon as you can.

Power Indicator Red: 0%-10% power left, please charge immediately.

While charging, the power indicator is blinking.

Power indicator will stop flash and stay bright when fully charged.

**\*CAUTIONS:**

1. Do not immerse the charging base in water or any liquid.
2. Do not use or place charging base in a high temperature environment.
3. Recommend use of power adapter with CE, UL or local safety standard 5V DC 1A output (need to purchase by the user).
4. Charging protection design adopted. Do not charge on the base for too long in order not to hurt the battery use life.
5. Do not use while charging.

**e. Cleaning**

Clean the brush head and keep it ventilated.

Whole body washable  
Able to wash under running water.

Change brush head every 2-3 months.

### III. Specifications

|                       |                  |                                               |
|-----------------------|------------------|-----------------------------------------------|
| Electric Brush Handle | Product Model    | DR.BEI C1                                     |
|                       | Rated Voltage    | 3.7V DC                                       |
|                       | Rated Power      | 2W                                            |
|                       | Waterproof Level | IPX7                                          |
|                       | Charging Time    | At least 4 hours from empty to fully charged. |
| Charging Base         | Rated Input      | 5V DC 1A                                      |
| Date of Production    |                  | Printed on the package.                       |
| Place of Production   |                  | China                                         |

## IV. Troubleshooting

| Problems                                                                           | Reasons                                                                                                                 | Solutions                                                                                                                                          |
|------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------|
| Brush head cannot attach to power shaft correctly.                                 | Stuff clog between the brush head and power shaft;<br>Brush head is not genuine or not matched;<br>Incorrect orienting. | Remove the blocking stuff;<br>Purchase genuine brush head recommended by manufacturer;<br>Align the parts correctly according to the instructions. |
| Cannot power on.                                                                   | Low power;<br>Product in charging.                                                                                      | Charge the product;<br>Remove brush from the charging base.                                                                                        |
| Weak brush head vibration.                                                         | Incorrect or loose attachment of brush head;<br>Stuff clog between the brush head and power shaft.                      | Attach correctly;<br>Remove the blocking stuff.                                                                                                    |
| Handle cannot attach to charging base.                                             | Stuff stuck at the bottom of handle;                                                                                    | Remove the blocking stuff.                                                                                                                         |
| The brush handle not showing charging status while attaching to the charging base. | The charging base failed to connect to the power supply.                                                                | Check the power supply.                                                                                                                            |
| Power indicator does not work.                                                     | Low power or used up.                                                                                                   | Charging correctly will soon restore it correctly.                                                                                                 |

## V. After-sales Service

If malfunction in DR.BEI Sonic Electronic Toothbrush Performance Failure List occurs, contact local authorized distributor or retailers.

### DR.BEI Sonic Electronic Toothbrush Performance Failure List

| Name                                     | Performance Failure                                                       |
|------------------------------------------|---------------------------------------------------------------------------|
| DR.BEI<br>Sonic Electronic<br>Toothbrush | Cannot turn on or off.                                                    |
|                                          | Brush head cannot connect to power shaft when there is no foreign matter. |
|                                          | Abnormal vibration of brush head.                                         |
|                                          | Power indicator does not work.                                            |
|                                          | Cannot charge the appliance.                                              |

## WARRANTY EXCLUSIONS

1. Damage caused by unauthorized repair, misuse, collision, neglect, abuse, liquid intake, accidents, alterations or incorrect use of product accessories that is not recommended by the manufacturer, or tear and alter labels or anti-counterfeiting marks.
2. No longer in the validity period of warranty service (repair, refund, replacement).
3. Damage caused by force majeure.
4. Failure to comply with the malfunction listed in DR.BEI Sonic Electronic Toothbrush Performance Failure List.
5. Failure of the product and its accessories in the DR.BEI Sonic Electronic Toothbrush Performance Failure List but caused by users.

## VI. Environmental Protection

When DR.BEI Sonic Electric Toothbrush reaches its useful life, you may need to remove the internal batteries and recycle them properly according to national or local laws and regulations. If you are unable to remove the battery, we strongly advise you to take your product to an official collection point designated by government.

EU Declaration of Conformity

We Wuxi Qinghe Xiaobei Technology Co., Ltd., hereby, declares that this equipment is in compliance with the applicable Directives and European Norms, and amendments.

### WEEE Disposal and Recycling

Correct Disposal of this product. This marking indicates that this product should not be disposal with other household wastes throughout the EU.

To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.

To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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方網站：[www.doctorbei.com](http://www.doctorbei.com)

Manufacturer: Wuxi Qinghe Xiaobei Technology Co., Ltd.

Address: Room 1106, 11th Floor, NO.2 of Building 530,  
Binhu District, Wuxi City, Jiangsu, P.R. China

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Адрес: Room 1106, 11th Floor, NO.2 of Building 530,  
Binhu District, Wuxi City, Jiangsu, Китай

почтовый ящик: [service@doctorbei.com](mailto:service@doctorbei.com)

メーカー：Wuxi Qinghe Xiaobei Technology Co., Ltd.

アドレス：Room1106, 11thFloor, NO. 2 of Building 530,  
BinhuDistrict, WuxiCity, Jiangsu, P. R. China

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제조업체: Wuxi Qinghe Xiaobei Technology Co., Ltd.

주소: Room 1106, 11th Floor, NO. 2 of Building 530,  
Binhu District, Wuxi City, Jiangsu, P.R. China

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For more information please go to DR.BEI official website:

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[www.doctorbei.com](http://www.doctorbei.com)

Scan to connect and  
discover more on Facebook.