artel

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

Floor scales

ART-EB-2802

PICTURE 1: TOP VIEW OF PRODUCT (FRONT SIDE)

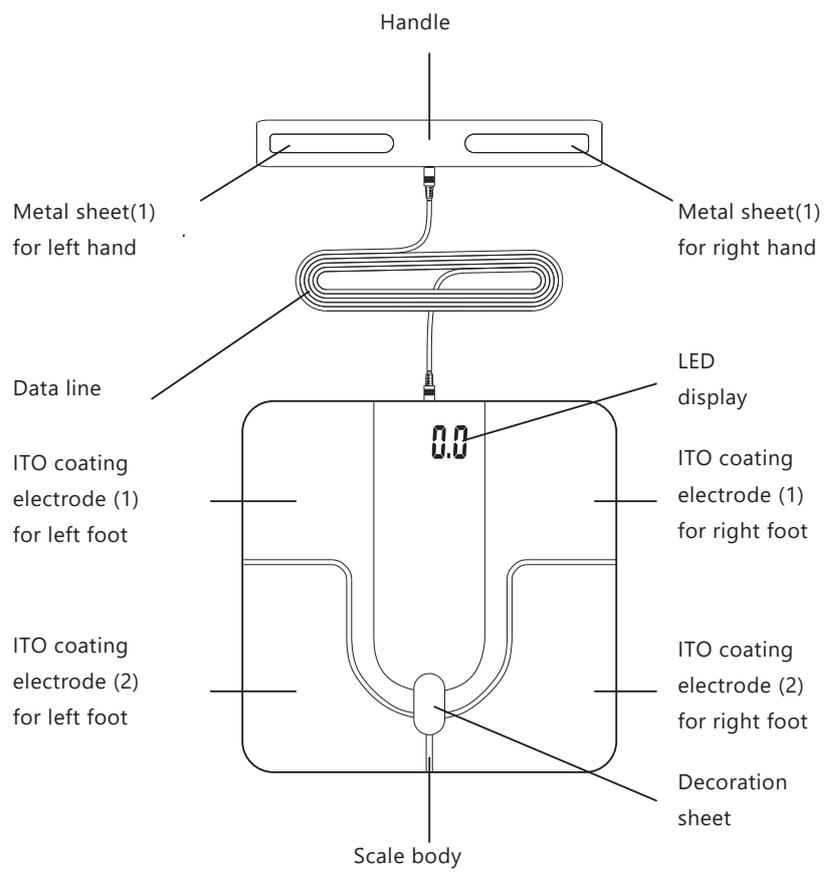PICTURE 2: TOP VIEW OF PRODUCT (BACK SIDE)

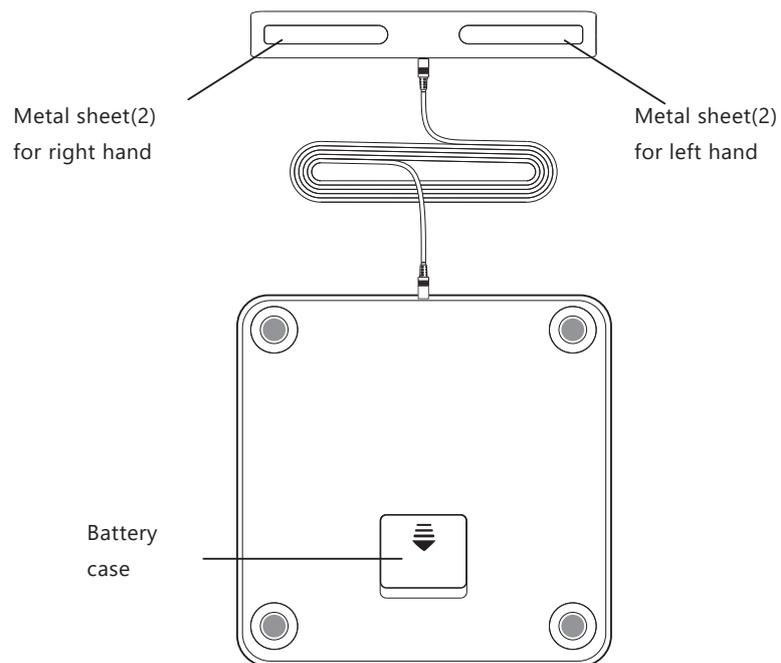
1. BATTERY INSTALLATION AND CHANGING

4x1.5V AAA batteries, please install as shown below:

First Use: Open the battery cover on the back of scale and take out the insulating strip.

Change battery: Open the battery cover on the back of scale, then install the batteries properly according to the pictures shown below.

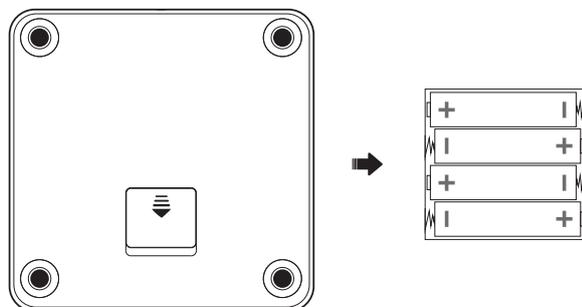Note: After battery installation, please close the battery cover. Please use corresponding battery. Remove battery if it is not used for a long time (more than three months).

2. DOWNLOAD AND SET UP THE APPLICATION

Search "MovingLife" in "App Store" or "Google Play" to download the Application. Or in another most convenient way by scanning the QR code below.


Tips: APP is kept updating, so it may differ from showed on manual. For the final and accurate information, please check on App.

3. REGISTERING AND CONNECTING


Open APP "MovingLife", connect the device with APP and submit your personal profile. (Before connection, please make sure Bluetooth is turned on and scale LCD is lighted up.)

Tips: Please make sure your setted personal information is correct, in order to get accurate measurements.

The LCD is lighted up

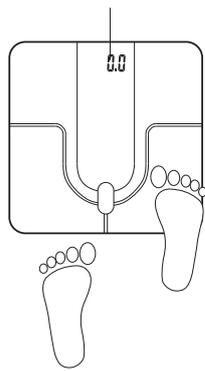Step gently to turn on the scale

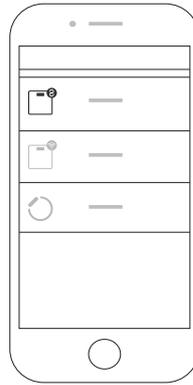
Connect the scale with APP

4. HANDLE INSTALLING AND MEASURING

• Installing


Connect the scale body and handle with data line as shown in picture.

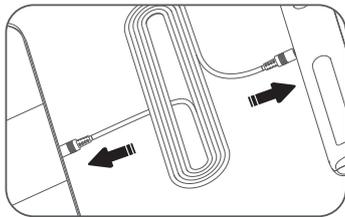
• Body composition measuring and analyzing

1. Place the scale on hard flat floor, open the APP on phone.
2. Lift up the handle, stand slightly onto the 4 ITO electrodes on the scale with bare feet; hold the 4 metal sheets on the handle with both hands, follow the tips shown on the handle, adjust the height of both hands and keep your arms unbend. (keep your arms away from your waist.)

Hold the 4 metal sheets on the handle with both hands , try your best to reach out your arms unbend. (keep your arms away from your waist.)

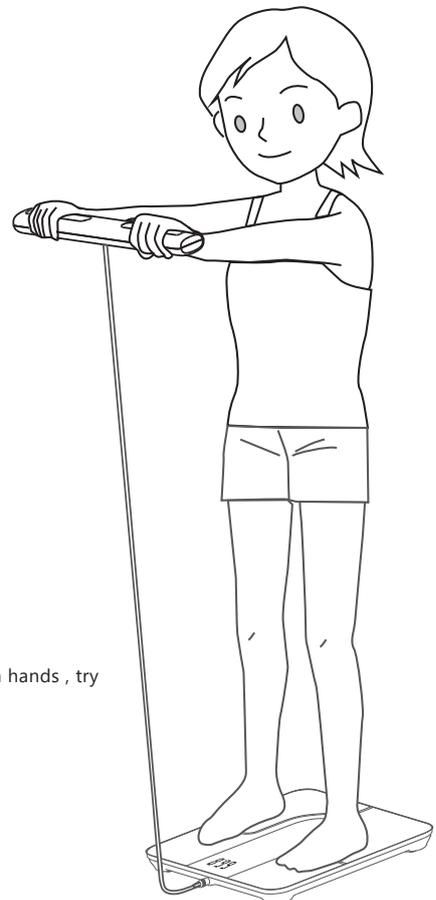- When the weighing result on the display is stable and flashes, the scale starts to analyze the body composition. After the progress bar is full, all the body composition data will be shown on APP and saved automatically. The display on the scale will only show weight, body fat and water rate circularly for 3 times, then the scale will be auto off.z

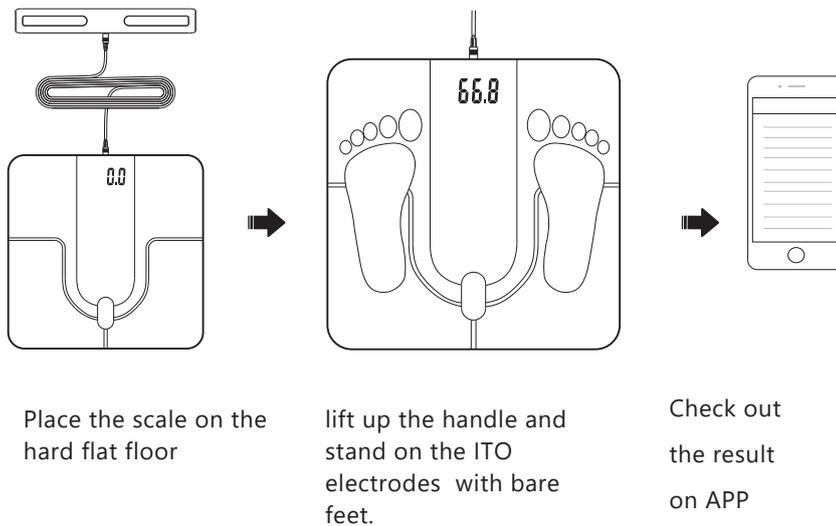
Tips:

- In order to ensure the measuring accuracy, please hold the 4 metal sheets on the handle tightly, stand on the scale with bare feet and keep both feet on the ITO electrodes. Please wait patiently until the measuring finished which takes about 20 seconds.
- If handle is not used, body composition can not be analyzed, only body fat and water rate can be measured.
- All the measuring should be finished with APP, if smart phone is not available, only weighing function is workable.

5. MULTIPLE USERS SUPPORTED

Multiple users blinding in one scale is supportable. Maximum 15 groups history data of each user can be stored. Date will be automatically -synchro-nizing to the mobile when APP is on.

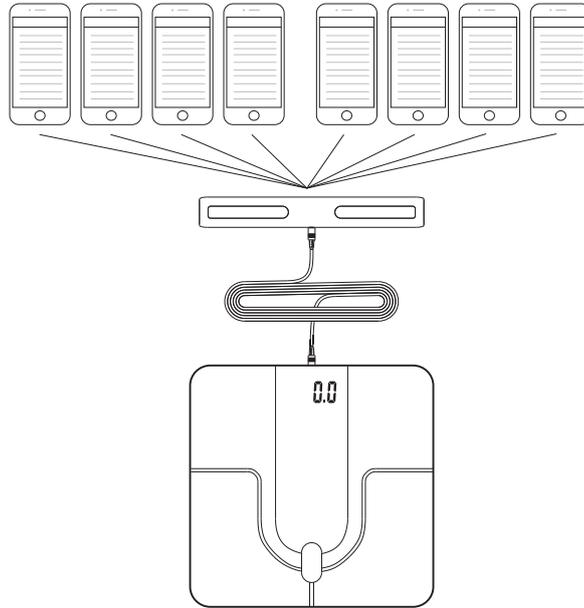Tips: Final user amount depends on the product.

6. DISPLAY ICONS

Measuring

Low battery
power

Over loaded

Re-measuring
needed


Body fat rate
is low

Body fat rate
is high

7. ABOUT THE UNIT

- To set the unit on APP is suggested.
- You can change the unit by the UNIT switch button on bottom of the scale (If the scale is with such button).
- The factory unit setting differs from different sales destination. Please refer to the actual display.

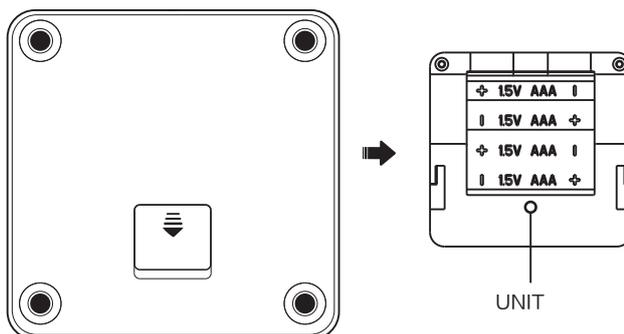BRIEF INTRODUCTION

This scale is a high precision body data analyzing instrument through weighing, 8 electrodes, 6 channels technology. It can transmit your weight data to smart phone via Bluetooth technology, as well as make graph of your weight results, from which you are able to keep track of the trend anytime.

Download the App on your phone to use bluetooth scale, you will start your first step of weight control.

ADVICE FOR USE & CARE

1. To ensure the accuracy of the data, it is recommended that each measurement be carried out at the same time of day. Don't use in strong electro-magnetic condition and field with intense shake .
2. If the scale leaves unused for long time, when first weighing, please tap on the scale to make LCD show zero and let the scale be off automatically. Then you can tap on again to start weighing.
3. To ensure the measuring accuracy, please take off the shoes and socks and other heavy belongs, keep your hands and feet clean, and then stand on the scale with bare feet and hold the handle with both hands.
4. Since the scale surface is smooth, please do not use the scale on slippery floor or stand on the edge of the scale; Always stand onto the central part of the scale and adjust your feet to touch the 4 ITO electrodes after you stand stably. Please keep standing firmly and balanced to prevent sliding.
5. Measuring results will be not accurate after strong exercises, overeating or in seriously dehydrated status.
6. Regarding the body composition analyzing, the measuring data may be not accurate for the below users:
 - Children below 18 or adults over 70
 - Professional sporters such as body builder
7. Below people are not suggested to use this product:
 - Pregnant women, Artificial dialysis patients, Bloated body status
 - People who are with heart pacemakers or other medical instruments inside the body are forbidden to use this product
8. When the batteries run out, please remove the batteries in time to prevent leakage of the battery lead scale body corrosion;
9. Do not drop, shock or strike the scale as it is a precise instrument.
10. Keep the scale clean. Do not use acid cleaning agent .

USER-FRIENDLY TIPS

The product can be used to weigh and calculate body fat % & hydration % after you set gender/ height/ age parameters ect. These measurements are only indicative for body fat status.

This product is not for medical treatment and can't be used in any medical institution. Growing children, pregnant women, those with symptom of edema and those on heart bypass treatment are not suggested to use the product.

FAQ

The software can't obtain weighing results from Smart scale.

- Please confirm the connection sign on the App interface to know the connection status.
- Confirm the Bluetooth is on in your device and the Smart scale is on within the valid range of signal coverage.
- Check whether the battery power is running low, please replace with a new battery.
- Relaunch the application.
- Turn off the Bluetooth and on again, then bridge the Bluetooth connection between scale and smart phone again.

Fail to connect to your device

- Please make sure only one device is setting up connection with the Smart scale.
- Please make sure your device is equipped with Bluetooth 4.0 or above.

LCD has no response

- Please first check the batteries and replace them if necessary
- Repairs may only be carried out by our authorized customer service or dealers.

The product picture is indicative only. Appearance specifications please in kind prevail.

artel

Designed by
Artel Engineering