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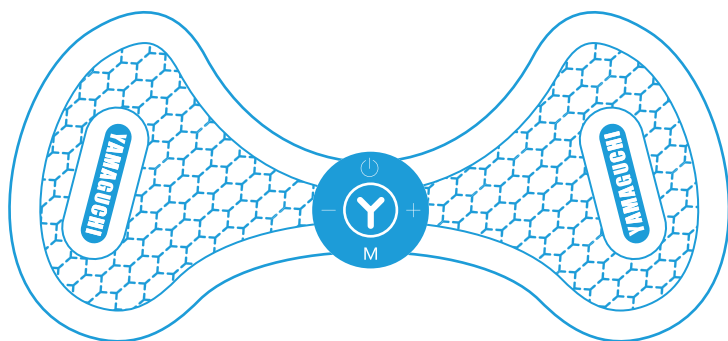
The **Yamaguchi NECK Trainer MIO** — is a next-generation portable electronic pulse massager with EMS technology developed to train men's and women's neck and lower back muscles. The device has 2 automatic modes with 15 adjustable levels of intensity to make trainings as much as possible comfortable and effective. Thanks to the unique and compact design, it's easy to use the device every day. Yamaguchi Neck Trainer MIO works from a powerful built-in battery and has a built-in 25 minute timer. The compact wireless design allows you to take it easy and to use it anywhere you wish. The Yamaguchi Neck Trainer MIO electronic pulse massager is an essential thing to improve blood circulation and restore lymph flow, relieve fatigue and muscle tension in the neck zone and lower back muscles, and helps to relax and quickly recover.

PRODUCT FEATURES

- Wireless use
- Suitable for the muscles of the neck zone and the muscles of the lower back
- Easy to operate
- 2 automatic modes
- 15 intensity levels
- Promotes relaxation and quick recovery
- It could relieve fatigue through promoting local blood circulation, relaxing local muscles and exercising the muscle tension
- Ergonomical design
- Electrode slice usage is 30-50 times
- Pocket size, easy to carry
- USB charge system
- Timer 25 min
- Ideal to use with a Yamaguchi wireless remote control

SAFETY INSTRUCTIONS

- Only use this device for the purpose for which it was developed and in the manner specified in the user manual instructions.
- This device is not designed to be used by persons (including children) with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless they have been instructed in the use of the device.
- It is forbidden to repair, modify, disassemble the device. If the faults exist, contact a certified service center.
- Never use the device while taking a bath or shower. If liquids penetrate the device, please avoid the further application. Contact to your certified service center.
- In case of diabetes or other diseases you should consult your doctor before applying the device.
- The use of device should be pleasant. If you feel pain or feel the procedure unpleasant, stop it.
- Ask your doctor for each kind of pain or disease, before using the device.
- Do not use the device while driving a car.
- If the device shows obvious damages at the main body or the electrodes, contact to your certified service center for repair.
- Do not use the unit near to equipment emitting a strong electric or magnetic field, such as radios, microwaves or

mobile phones. Otherwise this may impair the correct functioning of the unit.

- Do not use the unit near any heat sources, highly flammable substances, gases or explosives.
- If you wear jewelry or piercings make sure that you take it off before using the unit. Otherwise it may cause the possibility of injury.
- Prohibit using the device while sleeping.
- Do not use this device simultaneously with other devices which give electrical pulses to your body.
- It is recommended to use the device once a day. The session should not exceed 25 minutes.

Safety information about the electrode pad:

- Do not place the electrode slice to the damaged areas of the body.
- The device must not come in contact with water. Do not expose them to moisture.
- Take off the electrode pad before taking a shower / bath or swimming.

- Do not use this device directly after taking a shower or a bath.
- Make sure, that the skin parts are completely dry and clean before applying the electrodes.
- In case of skin irritation, stop using the device and consult a doctor.
- The contact of skin with electrodes must be very tight.
- The electrode pad can be used for 30-50 times.

In the course of continuous improvement of the device we reserve the right to make technical and design changes.


ATTENTION! This product is not intended to be used for medical purposes, including for the prevention, diagnosis, treatment, monitoring the state of the human body, medical research, rehabilitation, replacement, change of the anatomical structure or physiological functions of the body, preventing or interrupting pregnancy.

CONTRAINDICATIONS

It is not allowed to use the device in the following cases:

- if you have pacemaker or other internal medical electroapparatuses;
- hypersensitivity to the impulse current;
- pregnancy;
- menstruation;
- serious heart diseases;
- heart arrhythmia;
- urolithiasis disease, cholelithiasis disease and other GIT diseases;
- thrombophlebitis;
- exacerbation of chronic diseases;
- oncology;
- psychosocial disabilities;
- nephrotonia;
- purulent inflammation, open wounds, injury, burns on the skin;
- rosacea;
- varicosis;
- tuberculosis;
- Influenza, ARVI.

PRODUCT INFORMATION

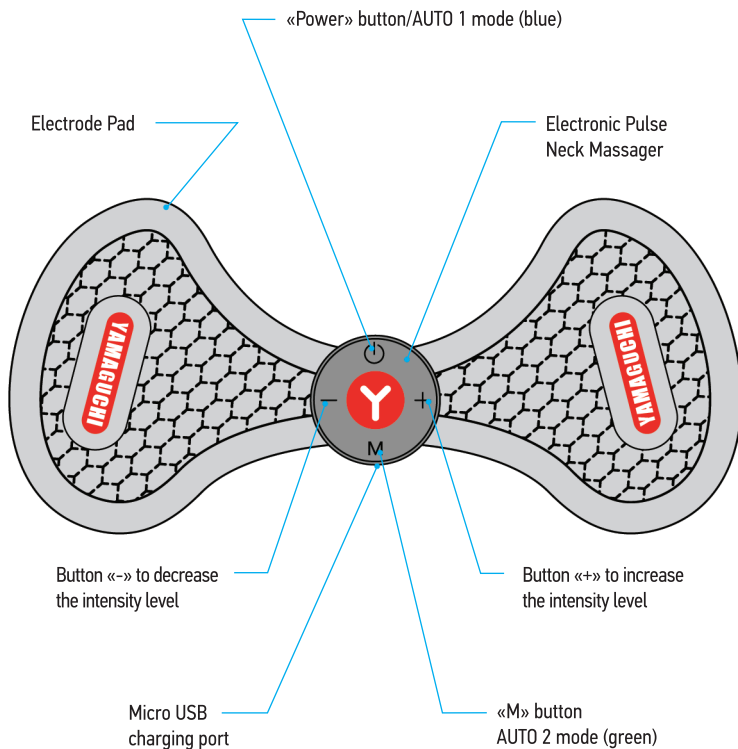PACKAGE CONTENTS

- Electronic Pulse Neck Massager x 1
- Electrode Pad x 1
- USB charging cable x 1
- User Manual x 1
- Package Box* x 1

*The package box is an integral part of the product appearance.

USER MANUAL

Read the user manual carefully before using this device, especially the safety instructions and keep the user manual for future use. If you give this device to another person, it is vital that you also pass on these manual for use.

To achieve the maximum effect it is recommended to apply before using scrubbing cosmetics to the skin, in order to make the electrode pad very tight.


Preparations before use.

Before using for the first time, charge the device using the USB charging cable that comes with the kit. Take the main body and attach the electrode pad to the main body, the device is ready for use.

Start to use.

1. Carefully remove the protective sticker from the electrode pad, without touching the gel coating by fingers, and stick the electrode pad to the muscles of your body.
2. Press the «Power» button and hold for 3 seconds to turn on the device. The blue indicator will light up and the AUTO 1 program will start up. Press the «+» button to increase or the «-» button to decrease the workout intensity level. There are 15 intensity levels. After 25 minutes the device turns off automatically.


AUTO 1 PROGRAM

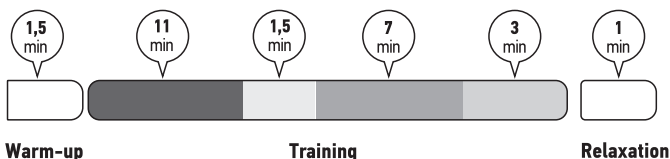3. Press the «M» button to turn on the AUTO 2 program, the green indicator light will turn on. Press the «+» button to increase or the «-» button to decrease the intensity level. There are 15 intensity levels. After 25 minutes the device turns off automatically.

AUTO 2 PROGRAM

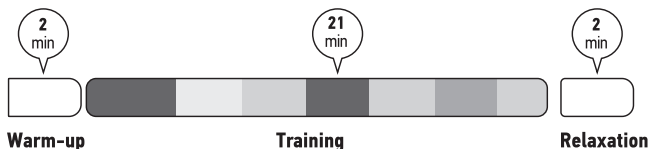4. With poor contact of skin with the electrode pad, the intensity level will automatically be reset to 0 level and after 30 seconds the device will turn off automatically. If the electrode pad is poorly attached to the skin, and has lost its stickiness, replace it with a new one. The use of one electrode pad is designed for 30-50 applications.
5. The device has a running time of 25 minutes by default. To turn off the device, press and hold the «Power» button for 3 seconds. Then gently remove the electrode pad from your body and stick the protective sticker for care and storage.

Charging the device.

Connect the USB charging cable one side to the power supply, the other to the charging port on the main body, the red indicator will light up on the display. When the device is fully charged, the blue indicator will light up. Please use a power supply with DC characteristics of 5V / 200mA and more. The full charge takes no more than 1 hour. The device does not work during charging process.

TROUBLESHOOTING

FAULT	CAUSE	SOLUTION
It cannot work normally	The battery runs low or out	Charge the battery
	The button is damaged	Contact the certified service center
The buttons cannot work normally	The battery runs low or out	Replace the battery
	The button is damaged	Contact the certified service center
The buzzer cannot beep for prompt	The buzzer is damaged	Contact the certified service center
No pulse output	The battery runs low or out	Charge the battery
	The electrode pad is damaged or polluted without any viscosity or has a bad contact with skin	Replace with a new electrode pad
	The skin is too greasy	Wipe skin with alcohol or wet towel
	This product is damaged	Contact the certified service center

FAULT	CAUSE	SOLUTION
Strong irritation with skin becoming red	Whether the massage lasts too long	It is most suitable to massage for 25 minutes
	Whether the massage intensity is too high	Decrease the intensity to a suitable level
	Whether the electrode pad doesn't contact close to skin	Contact the electrode pad close to your skin
	The skin is allergic to the electrode pad	Check if your skin is allergic. If slightly allergic, please change the service part of the electrode pad or shorten the treatment time. If seriously allergic, please seek for desensitization treatment
	Whether the electrode pad is dirty	Clean the electrode pad

CARE AND STORAGE

- Before each new procedure, the electrode pad must be disinfected for hygiene purposes, and to remove the remnants of sweat and fat (for better conductivity). The adhesive surfaces may only be cleaned by holding them under smoothly running, lukewarm tap water. Use your finger to slightly wipe the adhesive side to get rid off the remaining water. Do not use towels to clean or dry the adhesive surfaces in order to keep the gluing ability.
- Keep the adhesive surface of the electrode slice clean and avoid contact with dust, liquids, sharp objects etc. The re-usable protection foil should be put back onto the adhesive surfaces after each use.
- Clean the unit from time to time with a soft, dry cloth. Never use strong cleaning agents, very hot water, alcohol, naphtha, thinner or gasoline etc.
- Never immerse the device or any component in water. Be cautious not to get any moisture inside of the device.
- Do not expose the device to direct sunlight; protect it against dirt and moisture.
- Keep the device in the storage box when not in use.
- Store the device in a clean, dry location.

TECHNICAL SPECIFICATIONS

Product Name	Electronic Pulse Neck Massager
Model	Yamaguchi Neck Trainer MIO
Power Supply	DC 5V, 200mAh
Frequency	35-100Hz
Pulse width	20-400ms
Timer	25min
Device Dimensions	23,6 x 10,7 x 0,2cm
Package Dimensions	26,5 x 13,5 x 2,5cm
Net Weight	50g
Gross Weight	150g
Ambient environment	5-40°C
Relative humidity	≤65%
Atmospheric pressure	700-1060hPa

YAMAGUCHI®

by YAMAGUCHI ELECTRONICS SOLUTIONS, CO. LTD.

www.yamaguchi.ru

ERC