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FACE SKIN CARE DEVICE YAMAGUCHI EMS FACE LIFTING

Yamaguchi EMS Face Lifting is an innovative device with EFace Skin Care Device

MS microcurrent technology for lifting the skin of the face and décolleté, providing non-surgical lifting and the effect of facial contouring.

Japanese microcurrent technology EMS MicroCurrent consists in artificial stimulation of facial muscle contraction with electric impulses of a certain strength and frequency in order to activate metabolic processes in tissues, strengthen facial muscles, improve skin firmness and elasticity, tighten the face contour and eliminate mimic wrinkles. The impact of EMS pulses on the skin of the face and décolleté is a new effective procedure for rejuvenating and improving the quality of the skin.

Beauty Light Therapy has a beneficial effect on the skin and activates the rejuvenating functions of the skin. Blue light (wavelength $460 \pm 20\text{nm}$) soothes the skin, improves sebum secretion, has an anti-inflammatory and anti-edema effect. Red light (wavelength $630 \pm 20\text{nm}$) activates the regeneration of skin cells, stimulates the

YAMAGUCHI

production of collagen and elastin, actively combats the signs of skin aging. The thermal effect on the skin in combination with red light therapy increases the production of collagen fibers, which contributes to an increase in skin turgor.

Exposure to the skin with galvanic currents ION +\ - stimulates the improvement of the appearance and quality of the skin, the skin becomes smoother, younger and more radiant due to better absorption of cosmetic creams and serums.

The unique anti-aging vibration massage of tissues effectively reduces expression and superficial wrinkles, improves cell metabolism, eliminates puffiness and strengthens the muscles of the face and neck.

The unique shape of the platinum-coated working edge allows for a snug fit to the surface of the skin of the face, which contributes to the maximum efficiency of the procedure.

The ergonomic design and compact size of the Yamaguchi EMS Face Lifting will allow you to take the device on any trip and travel.

ENG

KEY FEATURES

- An innovative device with EMS microcurrent technology for face and décolleté skin tightening
- Non-surgical lifting and the effect of facial contouring
- Japanese microcurrent technology EMS MicroCurrent strengthens facial muscles, improves skin firmness and elasticity, tightens the face contour and eliminates expression lines.
- Beauty Light Therapy (14 LEDs) activates the regenerative and rejuvenating functions of the skin
 - ✓ Blue light (wavelength $460 \pm 20\text{nm}$) - soothes the skin, improves sebum secretion, has an anti-inflammatory and decongestant effect on the skin
 - ✓ Red light (wavelength $630 \pm 20\text{nm}$) - activates the regeneration of skin cells, stimulates the production of collagen and elastin, actively fights against the signs of skin aging
- Exposure to the skin with galvanic currents ION $+\backslash$ - improves the absorption of cosmetic creams and serums
- Anti-aging vibration tissue massage effectively reduces wrinkles, improves metabolism, eliminates puffiness and strengthens the muscles of the face and neck
- Heat exposure of the skin up to 45° increases the production

YAMAGUCHI

of collagen fibers, which contributes to skin rejuvenation, smoothing of wrinkles

- The unique shape of the working edge of stainless medical steel with platinum-plated coating allows a snug fit to the surface of the face
- Suitable for all skin types
- Compact size, ergonomic design
- Wireless use
- Timer 10 minutes

The positive effect of regular use of the device:


- the skin is 62% firmer;
- skin tone and elasticity are increased by 45%;
- the visibility of mimic wrinkles is reduced by 2 times.

Visually you will note:

- clearer and more toned oval of the face;
- elastic and firm skin;
- reduction of mimic wrinkles;
- improving complexion and evening out skin tone;
- fresh and radiant look.

DETAILS

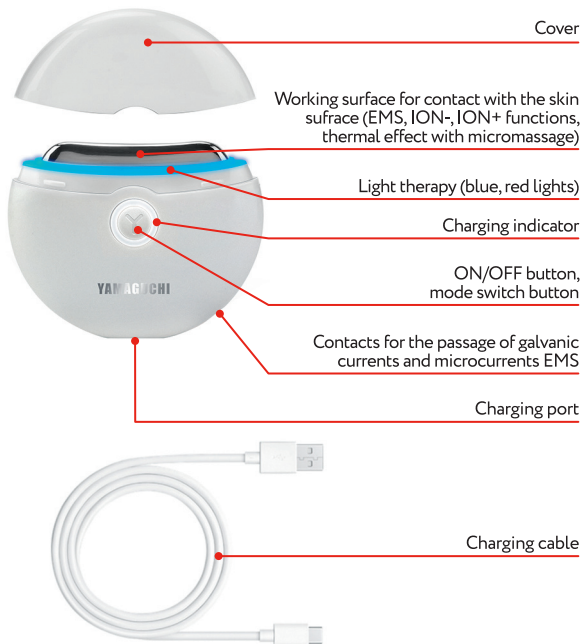YAMAGUCHI

USER MANUAL

Before using the device, carefully read the instructions for use, especially the safety instructions, and keep the instructions for use for future use.

Device charging

When using for the first time, fully charge the device.

When the battery level is low, the indicator will flash. The device needs to be charged.

- Connect the device to the adapter using the supplied USB cable. Use an adapter with DC 5V/1A or more. Connect the adapter to a 220V network, the charging process will begin. The charge indicator will flash white while charging. When the device is charged, the indicator will be on continuously.

Full charging takes no more than 4 hours.

Attention! If you do not use the lighting for a long time, the mirror should be charged every 3 months to extend the battery life.

Switching the device ON/OFF

1. To turn on the device, press and hold the ON/OFF button for 2 seconds.

2. The device has 3 operating modes:
 - 1 mode: light therapy with red light, exposure to galvanic currents ION + and thermal exposure to the skin up to 45 °
 - 2 mode: light therapy with blue light, exposure to galvanic currents ION-
 - 3 mode: red flashing light therapy, EMS microcurrent technology
3. Pressing the ON/OFF button twice activates the vibration micromassage mode. Pressing the ON/OFF button twice again deactivates the micromassage mode.
4. To turn off the device, press and hold the ON/OFF button for 2 seconds. The device will turn off automatically after 10 minutes.

User manual

1. Before using the device, completely clean your face and décolleté area.
2. Apply a small amount of any cosmetic product to your skin.
3. Press and hold the ON/OFF button for 2 seconds to turn on the device.
4. Select the desired mode by repeatedly pressing the ON/OFF button.
5. A total of 3 modes to choose from (in each mode you can turn

YAMAGUCHI

on additional micro-vibration by double pressing the ON/OFF button):

- 1 mode: light therapy with red light and exposure to galvanic currents ION +. Red light (wavelength $630 \pm 20\text{nm}$) activates the regeneration of skin cells, actively fights the aging process. Exposure to galvanic currents ION + in the cleansing mode deeply cleanse the skin at the cellular level. Thermal effect on the skin up to 45° improves skin microcirculation, opens pores, increases the production of collagen fibers;
- 2 mode: light therapy with blue light and exposure to galvanic currents ION-. Blue light (wavelength $460 \pm 20\text{nm}$) has an anti-inflammatory effect on the skin. Exposure to galvanic currents ION- actively nourishes the skin at the cellular level, promotes better absorption of creams and greatly enhances their effectiveness;
- 3 mode: flashing red light therapy and EMS therapy. Light therapy with red light (wavelength $630 \pm 20\text{nm}$) enhances blood circulation, helps to reduce wrinkles, makes the skin smoother and more elastic. EMS therapy acts on the muscles of the face with electrical impulses, strengthening them. Signals are sent to the nerve endings that control motor activity and force the muscles to contract in a natural way, thereby tightening the skin, the face acquires clear contours,

wrinkles are smoothed, and tissue "sagging" is eliminated. If during EMS facial therapy you feel a slight tingling sensation and muscle contraction, then the procedure is proceeding correctly;

When choosing third modes, it is necessary to apply a water-based gel or serum to the face for better current conductivity. We recommend using serum with low molecular weight hyaluronic acid Yamaguchi Hyaluronic Acid 2%. The serum improves skin elasticity, fills in already formed wrinkles, moisturizes the skin, accelerates regeneration, and eliminates minor inflammations.

6. Before using the device on your face, try it on your hand.

7. To activate the effect of galvanic currents and EMS-therapy,

YAMAGUCHI

it is necessary to touch the face skin with the working surface of the device, while touching the contacts on the sides of the device with your fingers.

8. Carry out, without pressing, the procedure along the massage lines. Do not linger on areas of the face.
9. **Upper third of the face.** Start moving from the eyebrow line to the hairline. This will help smooth out the crease lines on the forehead and make the skin look smoother and more even.

10. **Middle third of the face.** Sweep the device from the nose to the temples, grasping the area around the eyes with the edge of the working surface. Such movements reduce puffiness around the eyes and fine lines.

11. **Lower third of the face.** Swipe the device from bottom to top, starting from the chin to the zygomatic region. With this massage, the cheekbones will become more pronounced and beautiful contours will appear.

YAMAGUCHI

12. **Neck.** Massage from the décolleté to the earlobes. Such massage movements form a clearer oval of the face, reduce the severity of the "double chin", remove excess fluid and puffiness.

Attention! Do not run the metal work surface of the area too close to the eyes and lips, do not touch the mucous membrane of the lower eyelid and eyeballs.

13. After the procedure, press and hold the ON/OFF button for 2 seconds to turn off the device. The device will turn off automatically after 10 minutes.
14. After use, wipe the working surface of the device with a clean cotton pad moistened with chlorhexidine.

SAFETY INSTRUCTIONS

Please read this user manual carefully.

- Use the device only for its intended purpose, following the recommendations outlined in this user manual.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience and / or knowledge, unless they are under the constant supervision of those responsible for them or have received from these are instructions on how to use the appliance.
- Before first use, inspect the device for damage.
- Before using the device, consult a dermatologist and make sure that your skin condition allows you to use the device. During use, monitor the condition of your skin.
- Do not drop the device.
- If there is visible damage to the instrument case, contact a certified service center.
- It is forbidden to repair, modify or disassemble the device yourself. If there is a problem, contact a certified service center.
- If any abnormality occurs in the operation of the device (for example, if there is a smell of smoke, suspicious noise, etc.), immediately stop using it!

YAMAGUCHI

- It is not recommended to use the device on damaged skin, wounds and other sensitive areas of the face after skin peeling procedures. If discomfort occurs, discontinue use immediately and consult your healthcare professional.
- Do not use attachments not recommended by the manufacturer.
- After moving the device from cold to warm, it is not recommended to turn on the device for at least 1 hour.
- Do not expose the device to sudden temperature changes.
- Before using the device for the first time, fully charge it using the USB charging cable.
- Before the first use, test the operation of the device on the skin of the hand.
- Do not use the device on hypersensitive skin, consult a dermatologist.
- Do not immerse the device in water or rinse it with water.
- Use of the device should not cause pain.
- It is not recommended to use the device in the area of the eyes and lips.
- It is forbidden to use EMS-stimulation for people with implanted electronic and metal devices (with the exception of dental implants).

- It is not recommended to use the device with metal braces.
- The total time of the procedure should not exceed 10 minutes.
- After the procedure, wipe the working surface of the device with a damp or dry cloth.
- The device is intended for individual use.

Attention! In the course of continuous improvement of the device, we reserve the right to technical and design changes.

Note! This product is not intended for medical purposes, including prevention, diagnosis, treatment, monitoring of the human body, medical research, restoration, replacement, alteration of the anatomical structure or physiological functions of the body, prevention or termination of pregnancy.

YAMAGUCHI

CONTRAINDICATIONS

It is not recommended to use the device in the following cases:

- Purulent inflammation, wounds, abrasions on the skin;
- Presence of a pacemaker or other electronic and metal implants;
- Dermatitis, psoriasis;
- Patient feeling unwell (it is recommended to postpone the procedure);
- With the formation of boils, herpes, conjunctivitis on the skin;
- For eczema or rosacea;
- Immediately after chemical peelings, mechanical cleaning or laser procedures;
- If the skin of the face shows signs of an allergic reaction or it is irritated;
- Hypersensitivity to impulse current;
- Pregnancy and breastfeeding period;
- Oncological diseases.

STORAGE AND CARE

- After each use, wipe the device with a clean cotton pad moistened with chlorhexidine.
- Do not use harsh detergents to clean the device.
- Never use alcohol, thinner, or benzine.
- Do not leave the device in direct sunlight.
- Do not immerse the device in water or rinse it under running water.
- Store the device in a clean, dry place out of the reach of children.

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User manual × 1

Packing box* × 1

*The packing box is an integral part of the presentation of the product.

YAMAGUCHI

TECHNICAL SPECIFICATIONS

Name	Face Skin Care Device
Model	Yamaguchi EMS Face Lifting
Battery voltage	Li-ion 3.7V
Power	4W
Battery capacity	500mAh
Charging current	5V/1A
Electroplating function voltage	3V
EMS pulse frequency	2.5 Hz
Duration of EMS-impulses	195 ms
LED-light	blue ($460 \pm 20\text{nm}$), red ($630 \pm 20\text{nm}$)
Vibration frequency	8000 ± 2000 vibrations per minute
Heating	to 45°C
Timer	10 min
Charging time	4 h
Dimensions	$92 \times 92 \times 28$ mm
Weight	200 g

YAMAGUCHI

Japan

Product of "Yamaguchi Electronics Solutions"
Misumi, Nagato, Yamaguchi, Japan

www.yamaguchi.ru

EAC