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USER MANUAL

BODY MASSAGER
YAMAGUCHI JOINT CARE

CONTENT

Product information.....	19
Key features.....	20
Details.....	21
Control panel	22
Charging process	23
User manual.....	23
Safety Instructions.....	26
Contraindications.....	28
Troubleshooting.....	29
Storage and care	30
Technical specifications	31
Package contents.....	31

PRODUCT INFORMATION


Yamaguchi Joint Care is an innovative joint massager that combines unique functions to massage the knee and elbow joints of the human body. The massager has a complex effect on the joints, improving the elasticity of the ligaments and muscle performance, relieving pain in the joint area, accelerating the movement of lymph and blood flow, thereby improving the supply of tissues with nutrients and oxygen.

The device has 3 automatic massage programs, 3 levels of air-compression intensity levels, heating function, delicate vibration massage technology and has a timer for 10 minutes. The combination of air-compression massage with vibration function and a heating function will quickly relieve discomfort and pain in the joints, normalize blood circulation, quickly restore strength, achieve maximum relaxation of full body. Its compact, modern design allows you to use it with maximum health benefits wherever you want!

KEY FEATURES

- Innovative device for effective and gentle joint massage
- Suitable for knee and elbow massage
- The complex effect on the joints improves the elasticity of the ligaments and muscle performance, relieves pain in the joint area
- Joint Care Technology: 16 built-in 2D Push Air bags gently massage the joint and accelerate the movement of lymph and blood flow, improving the supply of tissue nutrients and oxygen
- Delicate vibration massage Shake-shake in the joint area relieves discomfort and pain in the joints
- 3 levels of heating function: 40°C, 45°C, 50°C
- 3 automatic massage programs
- 3 levels of intensity of air-compression massage
- Smart Memory function: the device automatically remembers the last used auto mode and intensity level
- Wireless use for 60 minutes
- Intuitive touch control panel
- Ergonomic design, compact size
- Timer for 10 minutes

DETAILS

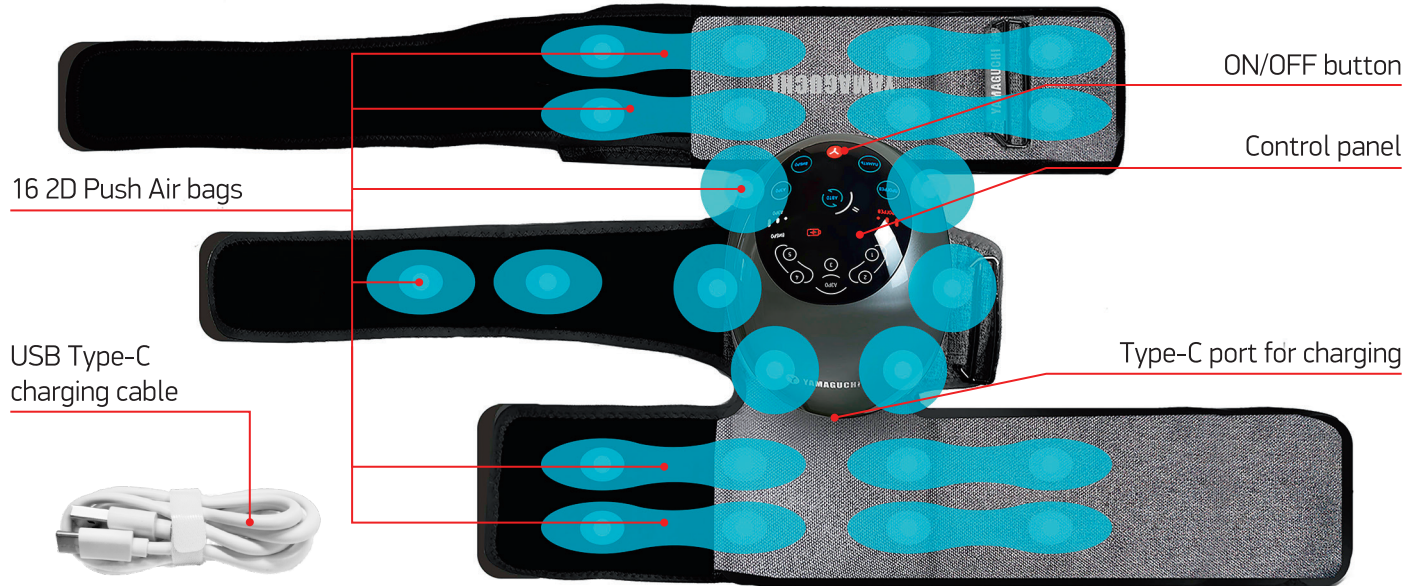CONTROL PANEL

Auto modes switching button

Smart Memory indicator

Heat function ON/OFF button/
1 of 3 levels of heating

Smart Memory ON/OFF button

ON/OFF Button

2D Push Air bags
working indicators

Battery indicator

Vibration massage indicator

Air-compression massage
ON/OFF button/
1 of 3 levels of air massage

Vibration massage
ON/OFF button


CHARGING PROCESS

In order to charge the device please connect the USB Type-C cable with one side to the power supply (not included), the other to the charging port, the red indicator will flash. When the device is charged, the indicator will stop flashing. Please use a power supply with characteristics DC5V 2A or more. It takes about 3 hours to fully charge the device. The massager does not work while charging.


*The average operating time of the device when the battery is fully charged is about 60 minutes.

USER MANUAL


1. Press the **Y** ON/OFF button located on the control panel in order to turn on the device.
2. Press the **AUTO** button repeatedly to select one of 3 automatic massage programs:
 - Auto mode I

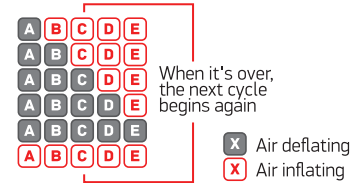- Auto mode II

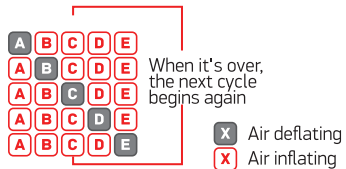- Auto mode III

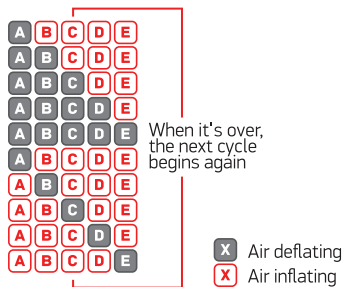3. Press the **AIR-COMPRESSION** Intensity level button repeatedly to select one of 3 levels:

- low
- middle
- strong

Note. Default setting: the LOW level.

Attention. We recommend turning on a low intensity level during the first massage session. With the help of textile fasteners, you can also fix a loose or snug fit, thereby regulating the strength of the massage.

4. Press the **HEATING** ON/OFF button to turn on/off the heating function during the massage.

Press the **HEATING** ON/OFF button repeatedly to select one of the 3 levels:

- level I - 40°C
- level II - 45°C
- level III - 50°C

Note. Default setting: the level I.

5. In order to remember the last using Auto massage program, the level of air-compression, heating and vibration, press the **MEMORY** button. The next time when you will use the massager, the last massage program will automatically start.

6. To enable/disable Vibration massage, press the **VIBRO** button.

7. To end the massage immediately, press the  ON/OFF button on the control panel. The massager will turn off and the air bags will gradually deflate.

SAFETY INSTRUCTIONS

Please read this user manual carefully.

- Only use this device for the purpose for which it was developed and in the manner specified in the user manual instructions.
- This device is not designed to be used by persons (including children) with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless they have been instructed in the use of the device.
- Before a massage procedure, please consult with your doctor and make sure that your health allows you to use the massager. During the massage, please take control of your health condition.
- The recommended massage time is 10 minutes per one session.
- Do not use the massager in the places with high humidity to avoid the electric shock.
- Do not use the massager close by any heat sources, highly flammable substances, gases or explosives.
- It is forbidden to use the device if the body of the device or its textile parts are wet.
- It is forbidden to repair, modify or disassemble the massager. If the faults exist, please contact a certified service center.
- Please stop use the massager immediately in case of any irregularity in the operation of the massager (for example, if there is a smoke or peculiar noise during use, etc.)!
- Do not cover the massager with a blanket during operation, as this may lead to overheating of the device.

- If you have any discomfort or feeling unwell during the massage, stop using the massager immediately and consult your doctor.
- Do not use the other attachments which are not recommended by the manufacturer.
- Do not use the device while sleeping.
- Do not get up or walk while using the device.
- Do not push foreign objects into the folds or gaps of the appliance.
- After you move the massager from cold to warm environment, the device is not recommended to use for at least 1 hour.
- Do not drop the massager.

CONTRAINDICATIONS

Please, consult your doctor before using the massager, if you have one of the following conditions:

- if you have pacemaker or other internal medical electroapparatuses;
- if you have different fractures, dislocations, ruptures of leg and feet muscles, tendon injuries and other external injuries;
- venous varicose;
- tenosynovitis;
- in cases of osteoporosis, thrombosis, severe arterial aneurysm, acute venous tumor, dermatitis, skin infections; in acute diseases;
- at high pressure and body temperature above 38°C;
- serious heart diseases;
- in the presence of diabetes mellitus;
- if you have contraindications to manual massage;
- oncology;

- in the presence of serious bone disorders, spinal fractures, sprains, seizures and acute pain;
- pregnancy and menstruation;
- if you have the injuries of the internal organs;
- children up to 12 years during the musculoskeletal system formation.

Attention! This product is not intended to be used for medical purposes, including for the prevention, diagnosis, treatment, monitoring the state of the human body, medical research, rehabilitation, replacement, change of the anatomical structure or physiological functions of the body, preventing or interrupting pregnancy.

TROUBLESHOOTING

Problem	Possible cause	Solution
Not operating	The product being charged does not work	This product is designed for safety protection and cannot be used when charging; after charging is completed, unplug the charging plug and press the power button to use this product
	The product has just finished a 10-minute massage and turned off the product	If you want to massage again, please press the power button again
Strange or loud sound	Vibrating sounds during massage are created by the product design	Vibrating sounds are absolutely normal
The button does not work	It may be that the lithium battery is insufficient	Lithium battery can be recharged
Product heat	The heating function is on	Turn heating function off
	The continuous working time of the product is more than 60 minutes	Product can be turned off and allowed to cool down

STORAGE AND CARE

- Store the massager in a cool, dry, well-ventilated place.
- Avoid touching the device with sharp or cutting objects, which may damage the tissue of the massager.
- Always turn off the massager before cleaning.
- Clean the body of the massager with a soft, slightly damp cloth.
- When cleaning the device, do not immerse the massager in water.
- Never use strong cleaning agents, very hot water, alcohol, naphtha, thinner or gasoline etc.

In case of malfunction of the device, immediately stop using it and contact YAMAGUCHI service center for testing, change and repair.

For safety reasons, all repairs of the device should be performed only by Yamaguchi service centers. The addresses and phone numbers of service centers can be found by tel. 8 800 333-12-81.

TECHNICAL SPECIFICATIONS

Product name	Body massager
Model	Yamaguchi Joint Care
Power	10W
Power supply	DC 5V/2A
Battery	li-ion battery 3,65V/2600mAh
Joint Care technology	yes
Air Bags quantity	16 pcs of 2D Push Air bags
Air Pressure	max 54 kPa
Shake-shake technology	yes
Vibration frequency	8000 rpm
Heating	40°C, 45°C, 50°C
Dimensions	18×11×14 cm
Net weight	1 kg
Timer	10 min
Temperature range	+ 5 °C ~ + 40 °C
Relative humidity	≤90%

PACKAGE CONTENTS

Body massager × 1
USB Type-C charging cable × 1
User manual × 1
Package box* × 1

* The package box is an important part of a product appearance.

Note. In the course of continuous improvement of the device we reserve the right to make technical and design changes.

YAMAGUCHI®
by YAMAGUCHI ELECTRONICS SOLUTIONS, CO. LTD.

ERC