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CONTENTS

Yamaguchi LED Light Therapy Mask.....	27
Product features.....	28
Positive benefits	
from using the LED Light face mask.....	31
Research results	31
Recommendations for use	31
Safety instructions	32
Pacakage contents.....	33
Product information	34
Remote controller description.....	35
User manual.....	36
Contraindications	41
Technical specifications	42
Care and storage.....	43

YAMAGUCHI

COSMETICS

LED-LIGHT FACE MASK YAMAGUCHI LED LIGHT THERAPY MASK


Yamaguchi LED Light Therapy Mask is an innovative wireless LED silicone mask for comprehensive facial skin care at home. The device combines unique cosmetology technologies, which together give a stunning effect, helps fight excessive sebum, pigmentation, acne, and helps narrow pores. The mask helps to improve the absorption of cosmetics, tightens and strengthens the oval of the face, improves skin nutrition and elasticity, thus increasing the production of collagen and elastin by skin cells, the facial contour becomes clearer and more prominent, the cheeks and double chin are tightened, wrinkles are eliminated, the skin is rejuvenated and acquires well-groomed blooming appearance.


LED Light therapy allows light waves of different lengths to penetrate into the cells of the deep layers of the skin, light is able to activate regenerative and rejuvenating functions, as a result of which the skin quickly becomes elastic and young. Skin looks smoother, fresher and radiant from within.

PRODUCT FEATURES

- Hi-tech Professional LED Skin Care device for personal home care
- The mask is made of hypoallergenic soft medical grade silicone
- 2 modes: auto mode and manual mode
- LED Light therapy: 4 LED light waves of different wavelengths for maximum results in facial rejuvenation:
 - **Red light** (wavelength 620 ± 10 nm) stimulates cells to renew, activates tissue regeneration. Suitable for flabby, atonic and tired skin with age-related changes. It activates the production of its own collagen and stimulates cellular metabolism, restoring complexion and toning the skin, helps fight wrinkles and scars. Red light penetrates to a depth of 6-7 mm, being absorbed in the epidermis and the skin itself (dermis), part of the light radiation reaches the subcutaneous fatty tissue.
 - **Blue light** (wavelength 460 ± 10 nm) has a soothing and anti-inflammatory effect on the skin, killing bacteria and viruses that cause acne and inflammation. The blue light treatment renews the skin and heals wounds, prevents scarring, shrinks enlarged pores, boosts cell metabolism, and makes skin smooth and radiant. Can be used for oily and combination skin. The depth of penetration of blue radiation into tissues is 0.5-1 mm.
 - **Orange light** (wavelength 590 nm ± 10 nm) directly affects the cells of the dermis - fibroblasts, stimulates the production of structural components of the skin: collagen, elastin, hyaluronic acid, enhances the skin's own energy, helps moisturize, smooth wrinkles, fights pigmentation, which allows apply it to mature, withered skin and to prevent premature aging. The depth of penetration of light into tissues is 2-4 mm.
 - **Violet light** (wavelength 700 ± 10 nm) is suitable for skin with reduced tone. Makes skin more elastic. Perfectly improves

complexion and fights deep wrinkles and swelling. Stimulates lymphatic drainage, improves blood circulation and reduces swelling. The depth of penetration into tissues is up to 10 mm, light radiation reaches the subcutaneous fatty tissue, which contributes to the breakdown of fats and the effective fight against edema.

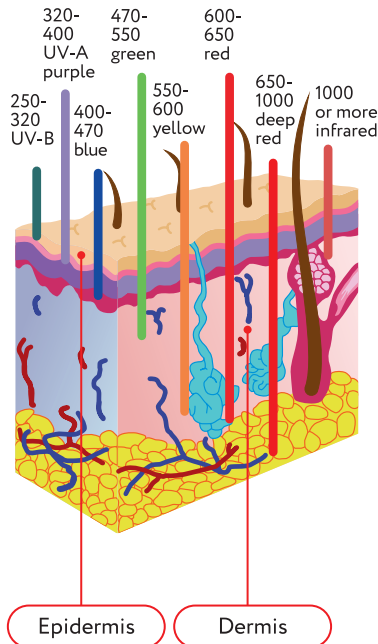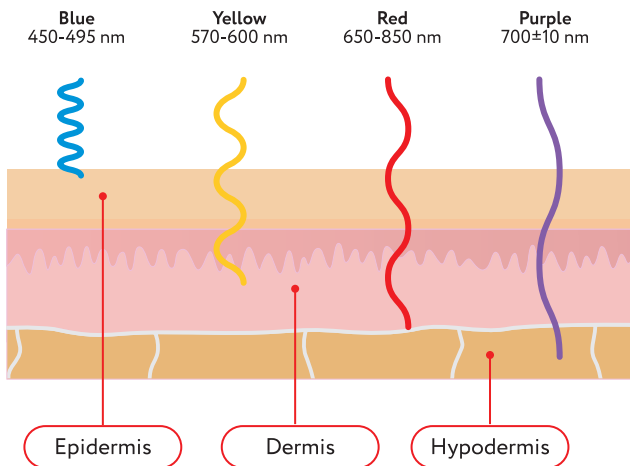- Convenient remote controller
- Additional head straps for a more comfortable fit of the mask
- Silicone protective rings around the eyes to prevent light waves from entering the eyes
- Silicone holes in the nose area, allowing you to breathe freely and wear a mask for a long time
- Easy to use
- Wireless use
- Ergonomical and compact design
- Simple Quick charge system
- Timer: auto mode 15 min, manual mode 10 min
- Suitable for all skin types and colors

POSITIVE BENEFITS FROM USING THE LED LIGHT FACE MASK:


- the skin becomes elastic, smooth and elastic;
- the oval of the face is tightened, the second chin disappears;
- stimulation of skin rejuvenation processes;
- elimination of pigmentation;
- smoothing wrinkles and nasolabial folds;
- the introduction of cosmetics into the deep layers of the skin, enhancing the effect of their impact;
- deep nutrition and hydration of the skin.

RESEARCH RESULTS:

- treatment with red LEDs reduces the number of wrinkles on the face by 38%;
- for 4-12 weeks, blue light therapy visually reduces the manifestation of acne and post-acne by 45-75%.

RECOMMENDATIONS FOR USE

- First week: Use once every 3 days, 15 minutes at a time.
- Weeks 2 and 3: Use once every 2 days, 15 minutes at a time.
- Fourth week: Use once a day, 15 minutes at a time.
- Then a break of 1-2 weeks.
- Then maintenance treatments 1-2 times a week, 15 minutes at a time.

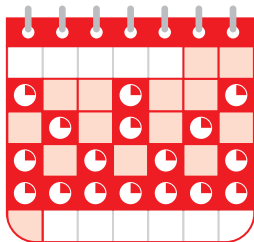SAFETY INSTRUCTIONS

Please read this user manual carefully.

- Only use this device for the purpose for which it was developed and in the manner specified in the user manual instructions.
- This device is not designed to be used by persons (including children) with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless they have been instructed in the use of the device.
- Before first use, inspect the mask for damage.
- Before use, consult a dermatologist and make sure that your skin condition allows you to use the device. During use, monitor the condition of your skin.
- Do not drop the device.
- In case of visible damage to the mask and cables, contact an authorized service center.
- Do not repair, modify or disassemble the device yourself. If there is a problem, contact an authorized service center.
- If any abnormality occurs in the operation of the device (for example, if you smell smoke, suspicious noise, etc.), stop using it immediately!
- It is not recommended to use the device on damaged skin, wounds, after skin peeling procedures. If discomfort occurs, discontinue use immediately and consult your physician.
- Do not use attachments not recommended by the manufacturer.
- After moving the device from cold to warm, it is not recommended to turn on the device for at least 1 hour.
- Do not expose the mask to sudden changes in temperature.
- Before using the led mask for the first time, fully charge it using the USB charging cable.

- Do not use the device for hypersensitive skin, consult a dermatologist.
- The instrument is not waterproof, do not store it in damp environments.
- Do not immerse the device in water or wash it with water.
- The device is intended for individual use.
- Before using the device, apply a water-based cream, gel or serum to the skin of the face and neck.
- In case of various diseases, it is necessary to consult a doctor before using the device.
- Do not use the mask while charging.
- Do not touch the mask, USB cable and USB adapter with wet hands while charging.


PACAKAGE CONTENTS

LED Light Face mask × 1
Head straps × 4
Remote controller × 1
USB-charging cable × 1
User manual × 1
Package box* × 1


* The package box is an integral part of the product appearance.

Note. In the course of continuous improvement of the device we reserve the right to make technical and design changes.

PRODUCT INFORMATION

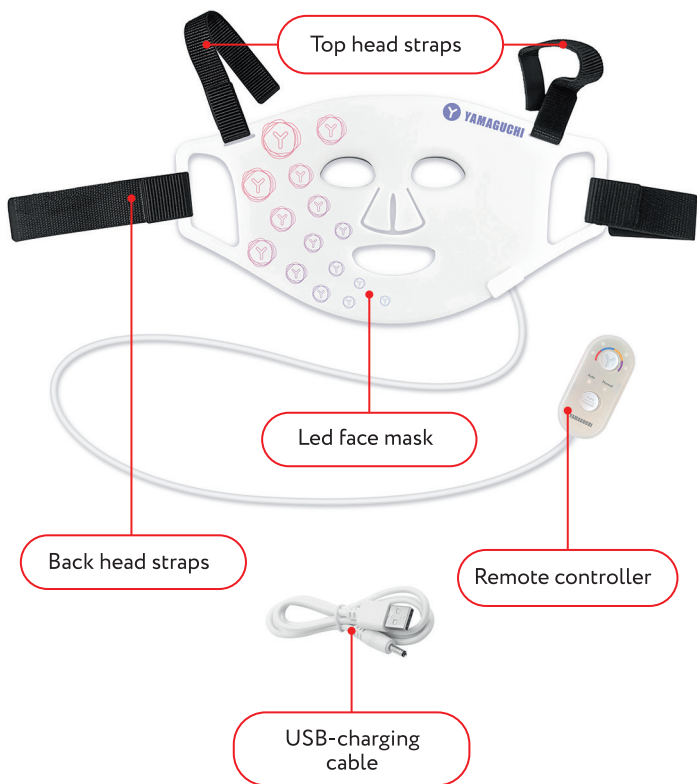REMOTE CONTROLLER DESCRIPTION

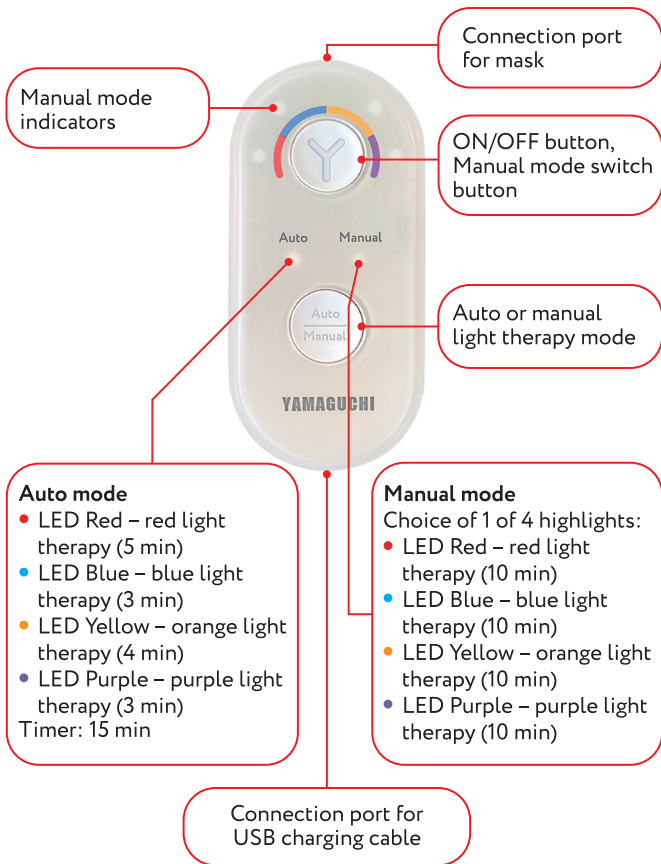
USER MANUAL

Read the user manual carefully before using this device, especially the safety instructions and keep the user manual for future use. Should you give this device to another person, it is vital that you also pass on these manual for use.

Charging process


1. When using for the first time, fully charge the device.
2. Connect the LED mask to the remote control to the appropriate port, connect the USB charging cable with one side to the remote control and the other side to the adapter (adapter not included). Use an adapter with DC 5V/1A specifications or more. Connect the adapter to a 220 V network, the charging process will begin.
3. During charging, the red light therapy indicator on the remote control will be lit. When the device is charged, all light therapy indicators will be on at the same time. A full charge takes about 2 hours.
4. It is forbidden to use the mask while charging.

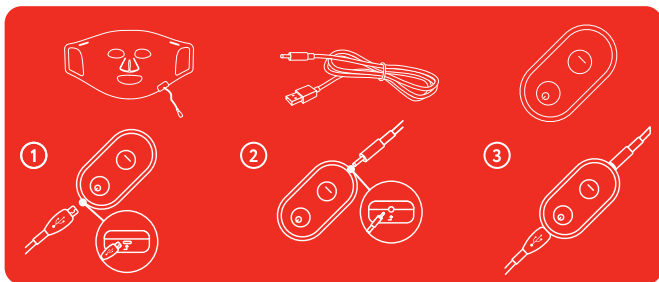Head straps

1. Included with the mask are special fixing straps for comfortable fastening of the LED mask on the head: 2 fixing straps for attaching from above and 2 fixing straps for fixing at the back.

2. Install the straps on the mask as shown in the figure.

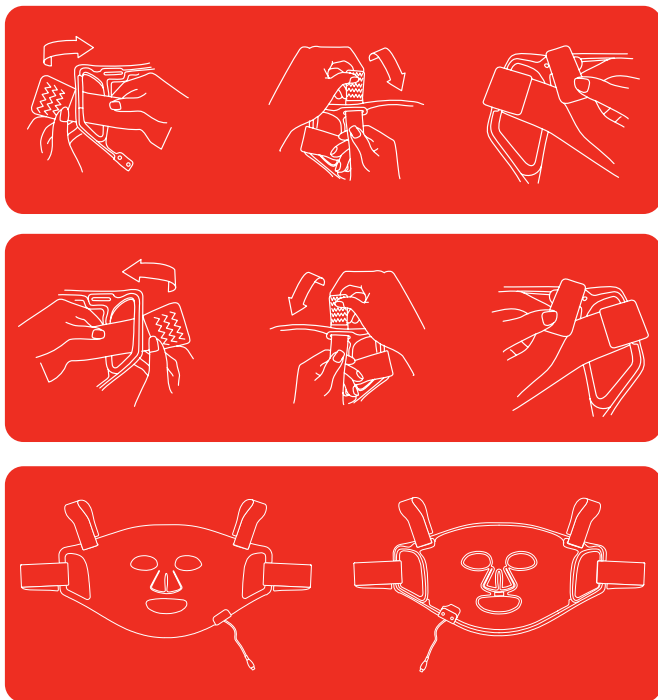Preparation for use

After installing the straps, connect the mask to the corresponding micro-USB connector on the control panel. The mask is ready to use.

User manual

1. Before using the mask, it is necessary to remove make-up of the skin of the face and neck, deeply cleanse the skin. Apply gel, serum or your favorite cream to the face and neck.
2. Put on the mask, fix the straps so that you are comfortable.
3. Turn on the mask by pressing the ON/OFF button. By default, the manual operation mode will turn on, the red LED will light up. Press the ON/OFF button again to select one of 4 light therapy options:

- **Red Light** (wavelength 620 ± 10 nm)
 - › Increases the synthesis of collagen and elastin;
 - › Accelerates metabolic processes;
 - › Improves the functional activity of cells;
 - › Corrects involutional skin changes;
 - › Helps reduce wrinkles;
 - › Reduces the depth of creases and folds;
 - › Increases skin density;
 - › Reduces the appearance of enlarged pores;
 - › Restores the level of skin hydration;
 - › Activates own defense mechanisms.

- **Blue light** (wavelength 460 ± 10 nm)
 - ▶ Soothes the skin;
 - ▶ Suppresses the development of inflammation;
 - ▶ Normalizes the work of the sebaceous glands;
 - ▶ Combats redness.

- **Orange light** (wavelength $590 \text{ nm} \pm 10 \text{ nm}$)
 - ▶ Stimulation of lymphatic drainage;
 - ▶ Elimination of puffiness and congestion;
 - ▶ Activation of microcirculation;
 - ▶ Antiseptic action;
 - ▶ Slowing down inflammatory processes;
 - ▶ Improvement of complexion;
 - ▶ Elimination of age spots;
 - ▶ Elimination of lack of oxygen.

- **Violet light** (wavelength $700 \pm 10 \text{ nm}$)
 - Reduced inflammation;
 - Antibacterial action;
 - Stimulation of regeneration;
 - Increased skin tone;
 - Reduced irritation, soothing effect;
 - Prevention of scarring in inflammation.

The mask will turn off automatically after 10 minutes
4. To enable auto operation, press the Auto/Manual button, the Auto indicator will light up. In this mode, all LED backlights are automatically changed:
- LED red light illumination (5 minutes)
 - LED blue light (3 minutes)
 - LED-illumination of orange light (4 minutes)
 - Purple LED illumination (3 minutes)
- The mask will turn off automatically after 15 minutes.
5. To turn off the mask, press and hold the ON/OFF button for 3 seconds.

ATTENTION! If you do not use the device for a long time, charge it every 3 months to extend the battery life.

CONTRAINDICATIONS

It is not allowed to use the device in the following cases:

- purulent inflammatory processes, wounds, abrasions on the skin;
- the presence of a pacemaker or other electronic implants;
- dermatitis, psoriasis;
- poor health of the patient (it is recommended to postpone the procedure);
- with the formation of boils, herpes, conjunctivitis on the skin;
- with eczema or rosacea;
- immediately after chemical peels, mechanical cleaning or laser procedures;
- if there are signs of an allergic reaction on the skin of the face or it is irritated;
- pregnancy;
- oncology;
- allergy to light;
- thyroid disease.

ATTENTION! This product is not intended to be used for medical purposes, including for the prevention, diagnosis, treatment, monitoring the state of the human body, medical research, rehabilitation, replacement, change of the anatomical structure or physiological functions of the body, preventing or interrupting pregnancy.

TECHNICAL SPECIFICATIONS

Product name	LED light face mask
Model	Yamaguchi LED Light Therapy Mask
Voltage	4,2 W
Battery capacity	1500 mAh
Battery voltage	3,7 V
Charging rate	5 V/ 1A
Light wavelength	Red light (wavelength 620±10nm) Blue light (wavelength 460 ±10 nm) Orange light (wavelength 590nm±10nm) Violet light (wavelength 700 ±10 nm)
LED diodes	100 pcs
Net weight	220 g
Dimensions	198×383×33,5 mm
Timer	Manual 10 min, Auto 15 min

CARE AND STORAGE

- Before each new procedure, it is necessary to disinfect the mask for hygiene purposes, as well as to remove sweat and grease residues.
- From time to time, wipe the mask with a soft dry cloth.
- Never use alcohol-based cleaners, thinner, or benzine.
- Do not submerge the mask or any of its accessories in water. Make sure that no moisture enters the device.
- Store the device in a clean and dry place.

In case of malfunction of the device, immediately stop using it and contact YAMAGUCHI service center for testing, change and repair.

For safety reasons, all repairs of the device should be performed only by Yamaguchi service centers. The addresses and phone numbers of service centers can be found by tel. 8 800 333-12-81.

YAMAGUCHI

Japan

YAMAGUCHI®

Product of "Yamaguchi Electronics Solutions"
Misumi, Nagato, Yamaguchi, Japan

EAC