{OOCOS | Aura

User Manual
One charge, One year, One toothbrush.

Contents

EN	1
DE	9
FR	18
ES	27
IT	36
PL	45
KR	54
JP	62
ZH	70
RU	78

Thank you for choosing Soocas oral care series: Soocas Aura Electric Toothbrush.

Please read this manual carefully and the correct method of usage for Soocas Aura Electric Toothbrush to fully understand its performance to better clean and care for your oral health.

Product Name: Soocas Sonic Electric Toothbrush

Model: PT1

Input Power: 5V==1A

Charging Time: 6 Hours

Charging Method: Type-C

Made in China

Manufactured by: Shenzhen Soocas Technology Co., Ltd.

Product Origin: Shenzhen, China

Address: Room 401, 402, 1101, 1102, 1103, 1104, Building 2 Chongwen Area, Nanshan iPark, No. 3370 Liuxian Avenue, Fuguang Community, Taoyuan Subdistrict, Nanshan District, Shenzhen, China

All products bearing this symbol are waste electrical and electronic equipment (WEEE as in directive 2012/19/EU) which should not be mixed with unsorted household waste. Instead, you should protect human health and the environment by handing over your waste equipment to a designated collection point for the recycling of waste electrical and electronic equipment, appointed by the government or local authorities.

Correct disposal and recycling will help prevent potential negative consequences to the environment and human health. Please contact installer or local authorities for more information about the location as the well as terms and conditions of such collection points.

Table of Content

Safety Tips	3
Product Introduction	4
Operating Guide	5
Charging Instructions	7
Environmental Protection	7
Troubleshooting	8

Safety Tips

When using this product for the first time, it is common to experience minor discomfort with some gum bleeding. This is because the gums are not quite used to the high-vibration cleaning method. It is recommended to use gentle mode for beginners and expects this discomfort to continue for 1-2 weeks upon initial usage. If bleeding persists for longer than 2 weeks, please stop using the product immediately and consult professional assistance from your dentist.


Warning

- For users with teeth correction, dental repair, dental implants, etc., please use this product under the guidance of a dentist or doctor.
- This product meets the safety standards of electromagnetic devices. For any issues related to pacemakers or other implantable devices, please consult a doctor or the manufacturer of the implant before use.
- This product is only designed for cleaning teeth, gums, and the tongue. Do not use it for any other purpose. If you feel any discomfort or pain, please stop using this product and seek professional assistance from a dentist or a doctor.
- Do not immerse the charging port in water or any liquid.
- Do not plug in the power supply if there is liquid on the charging port or power plug.
- This product has a built-in rechargeable battery. Do not throw it into the fire or charge it in a high-temperature environment.
- Do not knock or drop this product. If the product is damaged due to collision or falling, please stop using it.
- Do not modify, disassemble or repair this product by yourself.
- This device is not suitable for people with physical, sensory, or mental disabilities, and lack of experience and knowledge (including children) to use unless they are supervised or instructed by a person responsible for their safety.
- Children should not play with this product.
- Cleaning and maintenance should not be done by children without supervision.
- Use a safe and compliant power adapter for charging.

Cautious

- Please check for obvious damage before using the product. Do not use it if it is damaged or the machine has suffered a severe fall.
- Make sure that the charging port and charging cable are dry before charging; otherwise, it may cause an electric shock or short circuit.
- This product has a built-in battery, do not throw it into the fire.
- Do not place it in a high-temperature or damp place, and avoid direct sunlight.
- The charging base cannot be used in outdoor direct sunlight.
- Do not use the toothbrush when soaking in a bath or shower.
- Do not share the brush head with others. Do not place metal objects on the charging base while it is powered on.
- Note: The product and accessory illustrations in the manual are for reference only. Due to product updates and upgrades, the actual product may differ slightly from the illustration. Please refer to the actual product.
- This product is intended for home-use only. Do not use it in a commercial setting.

Product Introduction

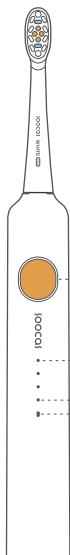
- Power Switch
- Mode Light
- Battery Reminder Lamp

Brush Head

Travel Case

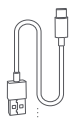
Charging Cable

Operating Guide


1. To Charge

- Please charge first before using it for the first time. Please use the 5V==1A charging adapter.

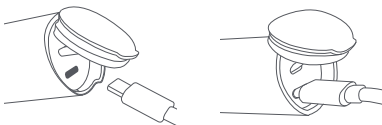2. Installing Brush Head

Ensure the bristle direction of the toothbrush head is aligned with the front of the toothbrush handle, then push the toothbrush head onto the motor shaft gently. Note that there will be a small gap between the toothbrush head and the handle to reduce friction between the two parts from the vibration of the motor.

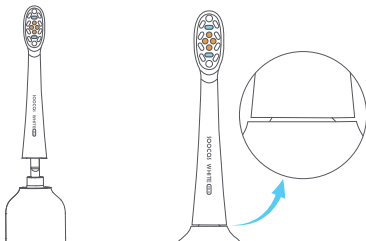3.To Turn On/Off and Pause

Key Function

1. Short press to start the device.
2. Short press the Power Switch within 15 seconds after startup to switch speeds. Short press the switch key within 15 seconds after startup to shut down.
3. In the startup state: Press and hold for 2 seconds to pause. If no action occurs within 30 seconds, the device will shut down, and if the Power Switch is pressed once again within 30 seconds, the device will once again start up and continue the countdown of the last action after the device is switched on.

Cleaning Modes

Soft: Gentle brushing, caring for new and sensitive users

Gum Protection: Comfortable and powerful, protects and improves gum health

Cleaning: Fulfills the daily requirements of most users

Brightening: Intensive cleaning power, restores the natural brightness

Polishing: High oscillation of the brush head for deep scrubbing of stubborn stains

Battery Reminder Lamp

Red light on: Low power, please charge in duetime

Red or green light flashing: Device is charging


Green light on: The equipment is fully charged

Zone Change Reminder

Every 30 seconds, it vibrates to pause and to remind replacing the brushing position.

4. Proper Brushing Techniques


1. Guide the brush head slowly from tooth to tooth with a 40-degree angle between the bristles and the gum lines.

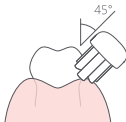2. Press the On/Off button for the preferred brushing mode.

3. The toothbrush will vibrate to indicate the user has to brush a new area of the mouth.

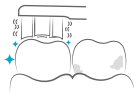4. Keep the brush head to the teeth closely without bending the bristles.

Charging Instructions

1. Soocas aura acoustic electric toothbrush has a built-in rechargeable battery, which takes about 6 hours to fully charge. The specific charging time may vary depending on the remaining power and charging environment.
2. Battery temperature should be 0 ~ 45 °C before charging.

Environmental Protection

When the product is at the end of its service life, according to local laws and regulations, one must remove the internal battery and undergo proper recycling and disposal procedures before the parts can be disposed of. If you are having trouble displacing the battery, please dispose of the product in accordance with the local waste clarification regulation.

Warning

- Do not throw batteries into the fire to prevent battery explosion or the release of toxic and harmful substances.
- Do not short-circuit the battery to prevent battery combustion and fire.
- Unless the intention is to discard the product, do not disassemble the product parts to prevent damaging the product.
- Strict safety procedures must be taken place. Before removing the battery, the power connection must be disconnected.
- Before discarding the product, the battery must be removed from the product.
- The product must be powered off when removing the battery.
- Batteries should be disposed of safely.

Remove the rechargeable battery

1. Use a flathead screwdriver to pry open the handle end cap.
2. Detach the cover and use a flathead screwdriver to pry open the inner motor. Pop out the motor.
3. Cut the red wires on the inner motor.

1. Insert a flathead screwdriver into the gap between the handle tail cover and the main body, and pry it upward until the tail cover opens for the left and right locking positions.

2. Press one end of the shaft against the desktop and apply downward force.


3. When performing the 2nd step, pry out the symmetrical two locking positions of the inner bracket with a flathead screwdriver, and apply downward force while holding the machine body to push out the movement.

4. Use diagonal pliers to cut it forcibly between the plastic and PCB plate at the point indicated by the arrow until the battery solder feet are cut.

5. Take out the battery and put it in the designated recycling place.

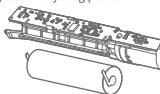Troubleshooting

Question	Probable cause	Solution
Unable to boot	Disabled when charging	Stop charging
Unable to boot, and red light flashes 3 times	Power is too low	Charge in time
Unable to boot, and red light flashes 6 times	Motor temperature is too high or too low	Please use it at a suitable ambient temperature (-20 °C ~ 60 °C)
Unable to install brush head	There are foreign bodies in the brush head hole	Remove foreign matter from the brush head hole
	There is a sticky foreign body on the rotating shaft	Remove sticky foreign matter
	Incorrect version of brush head	Please use the correct version of brush head
	The insertion direction is reversed	Please refer to the brush head installation instructions
Weak vibration of brush head	The brush head is not installed correctly	Reinstall the brush head

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