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Content

Product Features.....	14
Package Contents.....	14
Safety Instructions	15
Contraindications.....	16
Product Information.....	17
User Manual.....	18
Operating Modes Description.....	19
Troubleshooting	21
Care and Storage.....	21
Technical Specifications.....	22

Product Features

- Air compression foot massage
- Rolling foot massage
- 3 automatic massage programs
- Option of adjusting the intensity of air compression during automatic programs
- 2 types of rolling foot massage: reciprocating and rotating modes
- Suitable for home and office use
- Timer 15 min

Package Contents

- Foot Massager x 1
- Power Cord x 1
- Additional Removable Covers x 2
- User Manual x 1
- Package Box* x 1

*The package box is an important part of a product appearance.

Safety Instructions

Please, read this user manual carefully.

- Only use this device for the purpose for which it was developed and in the manner specified in the user manual instructions.
- This device is not designed to be used by persons (including children) with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless they have been instructed in the use of the device.
- Before the massage procedure, please consult with your doctor and make sure that your health allows you to use the foot massager. During the massage, please take control of your health condition.
- Before using the foot massager, please remove your shoes.
- Do not stand or seat on the foot massager and do not put the heavy objects on it.
- The recommended massage time is 15 minutes per one session.
- Do not move the foot massager during the massage procedure.
- Do not use the foot massager in the places with high humidity to avoid the electric shock.
- Do not use the foot massager close by any heat sources, highly flammable substances, gases or explosives.
- If the foot massager shows obvious damages at the body or power cord, please contact a certified service center.
- Make sure that the power cord is not twisted, tangled, bend and is not stretched in order to avoid damage or fire accidents.
- Do not touch the control panel of the foot massager with wet hands.
- It is forbidden to repair, modify or disassemble the foot massager. If the faults exist, please contact a certified service center.
- Please, stop use the foot massager immediately in case of any irregularity in the operation of the foot massager (for example, if there is a smoke or peculiar smell during use, etc.)!
- Do not pull the foot massager for the power cord.
- Do not use the massager on bare feet.
- Do not cover the foot massager with a blanket or a pillow when it is used, this can lead to device overheating.
- After each use, please turn off the power supply.
- During use, if you feel abnormal, please stop using the foot massager immediately and consult your doctor.
- Do not use the other attachments that are not recommended by the manufacturer.
- Do not insert the foreign objects into the folds or gaps of the foot massager.
- After you move the foot massager from cold to warm environment, the device is not recommended to be used for at least 1 hour.

Attention! In the course of continuous improvement of the device, we reserve the right to make technical and design changes.

Contraindications



It is not allowed to use the device in the following cases:

- If you have pacemaker or other internal medical electro apparatuses;
- If you have different fractures, dislocations, ruptures of leg and feet muscles, tendon injuries and other external injuries;
- Venous varicose;
- Osteoporosis;
- Influenza, ARVI;
- Hypertension;
- Serious heart diseases;
- If you have contraindications to manual massage;
- Oncology;
- Pregnancy and menstruation;
- If you have the injuries of the internal organs;
- Children up to 12 years during the musculoskeletal system formation.

Note.

Foot massager Yamaguchi Yume is not intended to be used for medical purposes, including for the prevention, diagnosis, treatment, monitoring the state of the human body, medical research, rehabilitation, replacement, change of the anatomical structure or physiological functions of the body, preventing or interrupting pregnancy.


Product Information

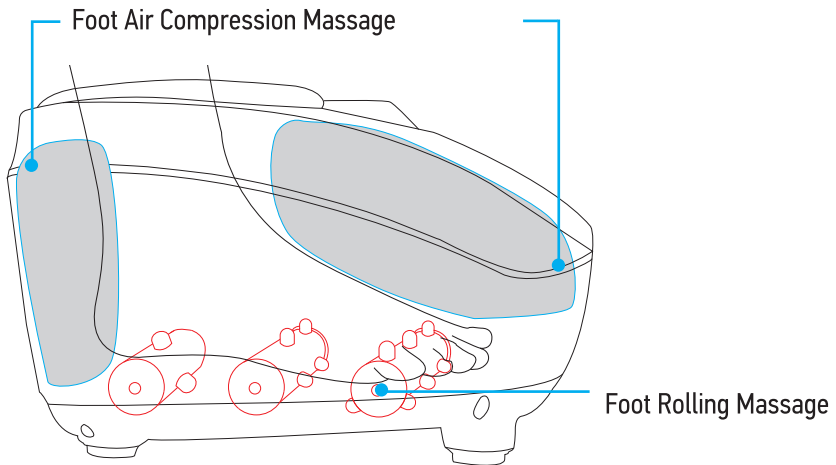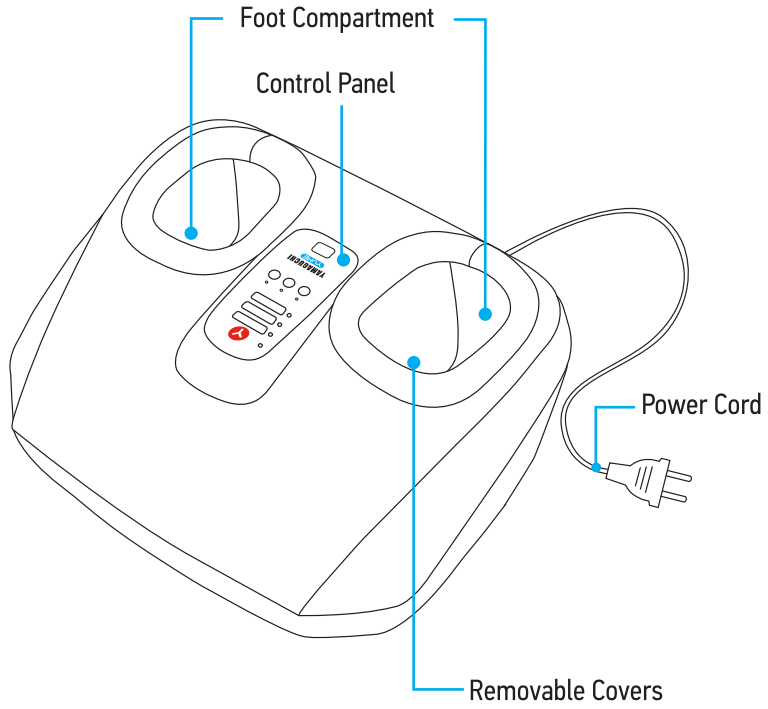User Manual


- Plug in Foot Massager to the power supply 220V.
- Make yourself comfortable in an armchair; put your feet in the holes for feet.
- Press **Y** the power button, located on the control panel, to turn on the device.
- To select one of the 3 automatic massage programs, choose the following modes:
 - DELICATE (Muscle relaxation, recovery their mobility and elasticity)
 - INTENSIVE (Prevention of heel spurs formation, water retention)
 - TONIC (Immune system reinforcement, acupressure)
- After 15 minutes of continuous operation the massager switches off automatically.
- Remove a plug from a socket.

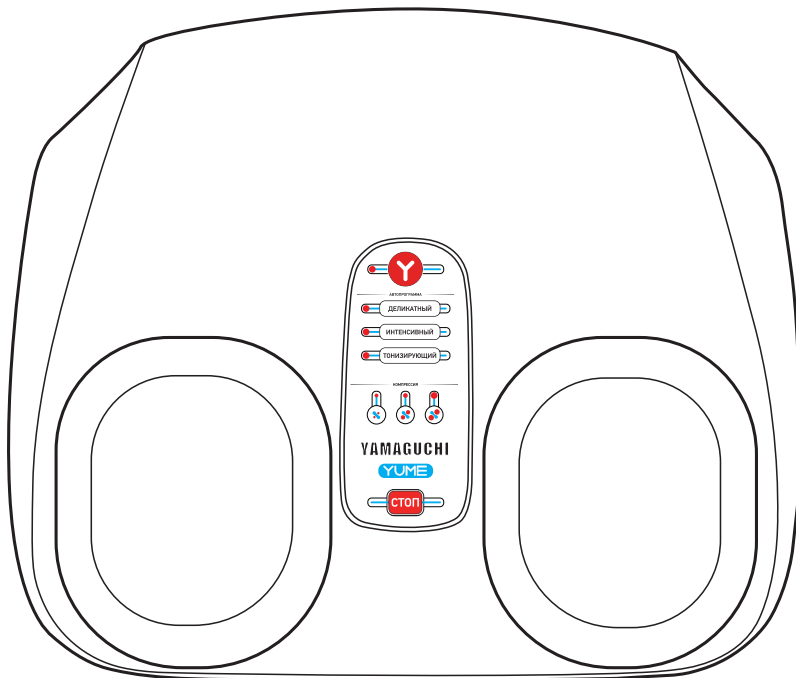Note. If you want to immediately stop the massage, press **СТОП** the button on the control panel of the foot massager.

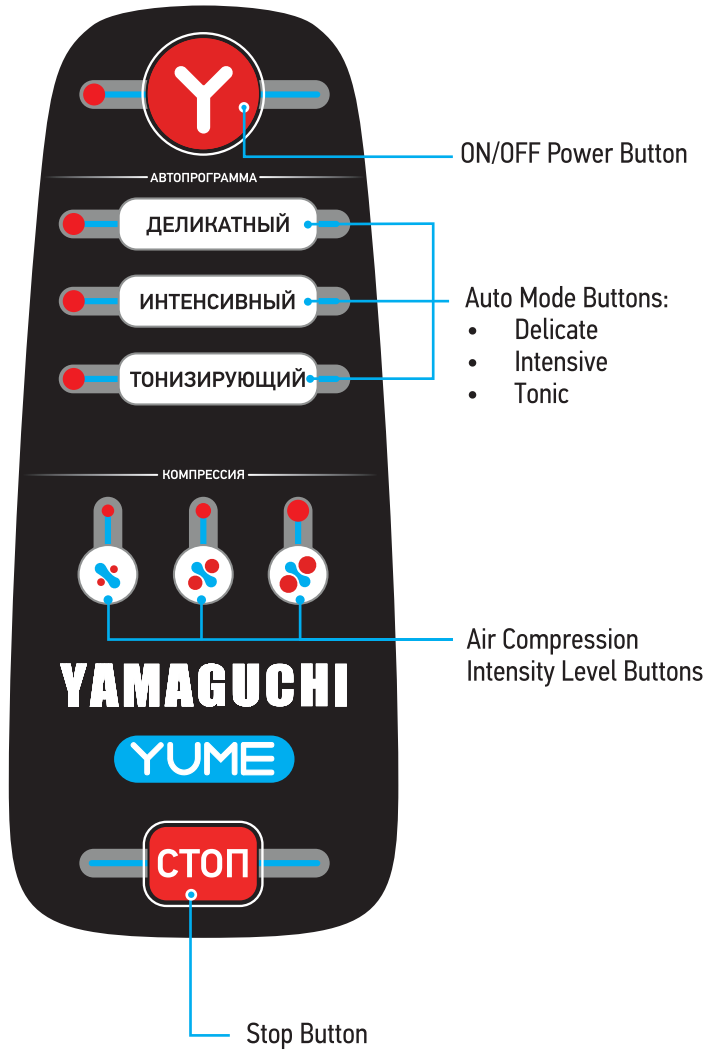Operating Modes Description

Control Panel

Massage programs	Description			
 <p>ДЕЛИКАТНЫЙ</p> <p>- Delicate -</p>	Press the  ON/OFF button to turn on foot massager then press Delicate button to select this automatic massage program. Delicate program has minimal rollers rotating speed:			
	Rollers Rotation Speed	Rolling Massage	Air Compression	Operating Mode
	Slow	+	+	Consecutive
 <p>ИНТЕНСИВНЫЙ</p> <p>- Intensive -</p>	Press the  “ON/OFF” button to turn on foot massager then press Intensive button to select this automatic massage program. Intensive program has maximum rollers rotating speed:			
	Rollers Rotation Speed	Rolling Massage	Air Compression	Operating Mode
	Fast	+	+	Mixed
 <p>ТОНИЗИРУЮЩИЙ</p> <p>- Tonic -</p>	Press the  “ON/OFF” button to turn on foot massager then press Tonic button to select this automatic massage program. Rollers rotating speed of this program is changing from minimal up to maximum:			
	Rollers Rotation Speed	Rolling Massage	Air Compression	Operating Mode
	Slow/Fast	+	+	Mixed

Press the “Air Compression Intensity” button to select the one of the 3 strength levels of air compression during the automatic programs:

LOW

MIDDLE

STRONG

Note. Default setting: the **LOW** level.

Attention! During the first procedure, we recommend to select **LOW** intensity level.

Troubleshooting

Problem	Cause	Solution
There is a clattering noise in the movement	This is perfectly normal caused by the working of rollers and air compression system	It is normal and no need to repair
The foot massager can't be started up	No power supply	Restore the electricity supply
	The power plug is loosened	Insert the plug into socket firmly
	The power switch is not turned on	Turn on the  power switch
Sudden STOP	Timer 15 minutes is off	Unplug the foot massager

Care and Storage

Care

- Before the cleaning, make sure the power supply is turned off and the plug has been pulled out from the socket.
- Use a soft cloth to wipe the surface of the foot massager and control panel.
- Never immerse the foot massager or any component in water. Be cautious not to get any moisture inside of the device.
- Never use strong cleaning agents, very hot water, alcohol, naphtha, thinner or gasoline etc.
- Removable covers can be washed in a washing machine at a temperature no higher than 30°C.

Storage

- Keep the foot massager in the original package.
- Store the foot massager in a dry and well-ventilated place.
- Avoid the contact of the foot massager with sharp or cutting objects that can damage the fabric removable covers.

Technical Specifications

Product Name	Foot Massager
Model	Yamaguchi Yume
Power Source	220 V / 50-60 Hz
Power Consumption	65 W
Dimensions (L x W x H)	48 x 40 x 23 cm
Net Weight	8,4 kg
Gross Weight	10 kg
Timer	15 min

YAMAGUCHI®
by YAMAGUCHI ELECTRONICS SOLUTIONS, CO. LTD