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CONTENTS

Head Massager Yamaguchi Headband.....	25
Key Features	26
The effect of using a massage headband.....	27
Details	28
Display	29
Massage mechanism Finger Touch.....	30
User manual.....	31
Charging	34
Package contents	34
Contraindications	35
Storage and care	35
Safety instructions	36
Troubleshooting	38
Technical specifications	39

HEAD MASSAGER YAMAGUCHI HEADBAND

Yamaguchi Headband is a wireless, compact and lightweight rim-shaped massager that combines unique point massage techniques, allowing you to cope with unpleasant sensations arising in the upper, occipital and frontal areas of the head. Vibration massage is carried out by the patented Japanese Finger Touch massage mechanism with 11 massage silicone fingers. Head massage improves blood circulation of the scalp, promotes oxygen saturation of the brain, increases brain activity and efficiency. Stimulating the hair follicles strengthens the hair roots and prevents their loss, which has a beneficial effect on hair growth. Massage effectively stimulates the acupuncture points of the head, thereby reducing fatigue, anxiety and tension.

The headband has 5 massage programs and 3 intensity levels, which separately and in combination, act as a measure of prevention of the main causes of pain localized in the head area. Additionally, the device can be used on other parts of the body — the neck and legs.

The functions are controlled via a remote control with an LED information display. Display is a removable battery as well, the remote can be charged with portable chargers autonomously.

The design of the product is made according to Japanese canons in carefully thought-out minimalism, which allows you to use the headband while travelling without additional transportation accessories. Compact size and light weight ensure convenient and easy use wherever you are comfortable.

KEY FEATURES

- Unique 2-in-1 device: elastic headband and head and neck massager
- The massager, in the form of an elastic band, can adapt to different sizes and features of the head shape, has many positions for wearing
- When adjusting the angle of wearing, the band covers 100% of the surface: the forehead, upper forehead area, parietal part of the head, back of the head, neck and legs
- The device copes with unpleasant sensations that occur in the upper, occipital and frontal parts of the head with unique acupressure techniques in the area of acupuncture points
- Patented Japanese Finger Touch massage mechanism: 11 independent silicone heads (“fingers”) perform unique massage combinations for quick relaxation and relief of fatigue
- The frequency of vibration 12,000 times per minute improves the blood circulation of the scalp and stimulates the acupuncture points of the head, improving brain activity and performance
- The high-frequency microvibrational Japanese motor operates in the range of up to 40 dB even at on the highest level of intensity, which will not cause discomfort to others
- 5 massage programs are designed to prevent the main causes of pain, stress and tension localized in the head and neck area:
 - ▶ **Wake up** — dynamic vibration massage (vibrations change every 0.15 seconds), will give you a boost of energy for the whole day.
 - ▶ **Relaxation** — the soft rhythm of massage combinations (vibrations change every 0.33 seconds) decreases psychoemotional tension and quickly calms down.
 - ▶ **Healthy sleep** — the slow rhythm of massage combinations and low intensity (vibrations change every 0.25 seconds) soothe, normalize the course of a night's rest and fight the symptoms of insomnia.
 - ▶ **Activation of hair growth** — unique acupressure techniques in the area of acupuncture points of the head improve blood circulation of the scalp,

stimulating the work of hair follicles and active hair growth

- ▶ **Neck massage** — 7 massage fingers in the middle of the rim perform a deep massage with high intensity, which allows you to relax the neck muscles.
- 3 levels of intensity with different massage programs increase the variability of the massage
- The control panel is equipped with a convenient LED display
- Japanese-style minimalistic design makes it easy to store and transport the product without additional accessories
- Joy Voice Assistant accompanies your commands
- Effective head and neck massage anytime, anywhere
- Wireless use
- Compact design
- Light weight: 120 g

THE EFFECT OF USING A MASSAGE HEADBAND:

- Activation of blood and lymph circulation;
- Improvement of the vascular tone of the scalp;
- Normalization of sleep;
- Stimulation of hair growth;
- Anti-stress therapy;
- Improves brain activity and performance;
- Increased concentration and attention;
- Immunity boosting.

DETAILS

DISPLAY

Charge indicator
Bluetooth
indicator

Massage mode
indicator

Intensity level:
low, medium, high

Wake up

Relaxation

Healthy
sleep

Activation of
hair growth

Neck
massage


MASSAGE MECHANISM FINGER TOUCH

It is a unique combination of massage techniques performed by 11 silicone massage fingers. The innovative mechanism creates the sensation of a professional massage with fingertips, interacting with delicate vibrations of different frequency and intensity on the acupuncture points of the head and neck.

The new massage technology allows a delicate action to eliminate tension and pain in the massaged area, normalizes blood circulation and improves metabolism, which allows you to enjoy a pleasant feeling of relaxation, recover and improve overall well-being quickly.

The flexible design provides a comfortable fit for any head and neck size.

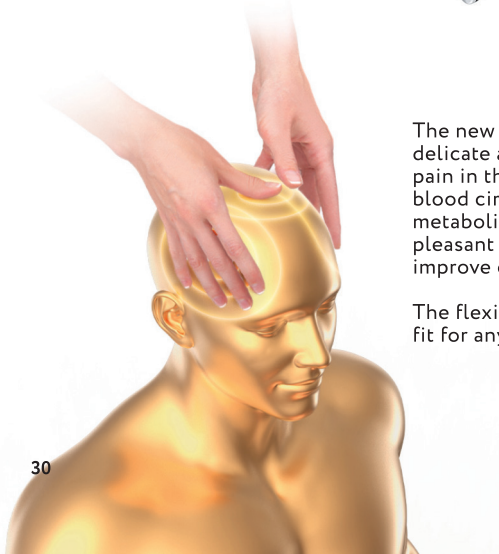USER MANUAL

1. Put the headband on top of the massaged area, ensuring that the silicone fingers fit the surface of the scalp or other parts of the body — neck or legs.

Forehead

Top of head


Middle of head

Back head

Neck

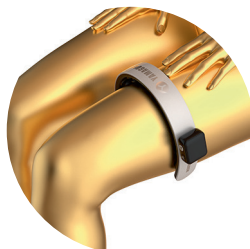Thigh

2. Press and hold the **Y ON/OFF** button to activate the massage band.
3. Press the **Y ON/OFF** button sequentially to select one of the 5 automatic massage programs. The current program indicator will appear on the information display.

	WAKE UP	Dynamic vibration massage (vibrations change every 0.15 seconds), will give you a boost of energy for the whole day
	RELAXATION	The soft rhythm of massage combinations (vibrations change every 0.33 seconds) decreases psychoemotional tension and calms down quickly
	HEALTHY SLEEP	The slow rhythm of massage combinations and low intensity (vibrations change every 0.25 seconds) soothe, normalize the course of a night's rest and fight the symptoms of insomnia
	ACTIVATION OF HAIR GROWTH	Unique acupressure techniques in the area of acupuncture points of the head improve blood circulation of the scalp, stimulating the work of hair follicles and active hair growth
	NECK MASSAGE	7 massage fingers in the middle of the rim perform a deep massage with high intensity, which allows you to relax the neck muscles

YAMAGUCHI

4. Consistently press the button to select the level of intensity of the massage effect. There are 3 levels available:
 - **Low intensity level** - delicate kneading, suitable for relieving tension
 - **Medium intensity level** - moderate-intensity shiatsu massage, relieves fatigue symptoms
 - **High level of intensity** - high-intensity acupuncture massage, helps against dizziness and headache
5. Press and hold the **Y ON/OFF** button to turn off the headband.
When the 15-minute timer ends, the device will turn off automatically.
6. **When is the best time to carry out the massage:**
 - ▶ it is recommended to massage the head daily, in the morning or in the evening;
 - ▶ during the day, when there is a need to relieve fatigue and pain;
 - ▶ before going to bed, to relax and relieve mental fatigue.

CHARGING

1. Please fully charge the massager before the first use. Charging cable (Type-C) is included in the package.
2. To charge, connect USB Type-C cable with one side to the power supply and the other side to the charging port. Use a power supply with the following characteristics, DC 5V/2A or more. The charging indicator will display the charging process on the display.
3. It takes 1.5 hours to fully charge the device. The average operating time of the device when the battery is fully charged is 520 minutes. When used once a day (for 15 minutes), the duration of use without recharging will be 35 days.

The device does not work during charging.

Note. If you don't use the device for a long time, charge it once a month to extend the battery life.

PACKAGE CONTENTS

Massage headband × 1
Remote control × 1
USB Type-C charging cable × 1
Storage Bag × 1
User manual × 1
Packing box* × 1

* The packaging box is an integral part of the product appearance.

CONTRAINDICATIONS

If you have any of the following conditions, consult your doctor before using the massager:

- electronic heart stimulator or other internal medical electrical devices;
- osteoporosis disease;
- acute diseases;
- hypertension;
- cardiovascular diseases;
- contraindications of manual massage;
- malignant tumors;
- pregnancy and menstruation;
- internal injuries;
- internal implants and prosthetics;
- age up to 12 years during the period of formation of the musculoskeletal system.

STORAGE AND CARE

- Store the massager in a cool, dry, well-ventilated place.
- The device is recommended to be stored at a temperature of +5 °C up to +40 °C and a relative humidity of not more than 90%. Protect the appliance from significant temperature changes.
- Avoid the device coming into contact with sharp or cutting objects that may cause damage to the case.
- Clean the body of the massager with a soft, slightly damp cloth.
- Do not immerse the massager in water while cleaning the product.
- Do not use abrasive cleaners, brushes, gasoline, kerosene, glass or furniture polishers, or paint solvents to clean the product.
- For long-term storage, it is recommended to charge the massager once a month.

SAFETY INSTRUCTIONS

Please read this manual carefully.

- Use the appliance only in accordance with its intended purpose, following the recommendations set out in this user's manual.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge, unless they are under the constant supervision of persons responsible for them or have received instructions from them on how to use the appliance.
- Before using it for the first time, inspect the massager and make sure that it is not damaged.
- Before starting the massage, consult your doctor and make sure that your health condition allows you to use the massager. During the massage, take care of your health.
- Do not use the massager in areas with high humidity to avoid electric shock.
- Do not use the massager near hot objects.
- Do not drop the massager.
- Never touch the massager with wet hands.
- Do not use the massager while driving.
- If there is visible damage, contact certified service center.
- Never repair, modify, or disassemble the device yourself. If there is a problem, contact the certified service center.
- It is not recommended to use the massager on joints, wounds and other sensitive parts of the body. If you feel uncomfortable or unwell during the massage, stop using the massager immediately and consult your doctor.
- Do not use any additional accessories that are not recommended by the manufacturer.
- After moving the massager from a cold room to a warm one, it is not recommended to turn on the device for at least 1 hour.
- It is forbidden to immerse the massager in water or other liquids while cleaning it.
- The recommended time for continuous massage is 15 minutes.

Note. In the process of continuous improvement of the device, we reserve the right to make technical and design changes that improve the operation of the product.

Attention! The Yamaguchi Headband is not intended for medical purposes, including diagnosis, treatment, monitoring of the human body, performing medical research, restoring, replacing, changing the anatomical structure or physiological functions of the body, preventing or terminating pregnancy.

If case of malfunction in the device, immediately stop using it and contact YAMAGUCHI service center for testing, change and repair.

For safety reasons, all repairs of the device should be performed only by YAMAGUCHI service centers. The addresses and phone numbers of service centers can be found by tel. 8 800 333-12-81.

TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
The massager doesn't work	Low battery power	Charge the massage band using the USB Type-C charging cable
	15-minute timer went OFF	Restart the timer by pressing the  ON/OFF button
	The battery is not installed in the battery compartment	Insert the charged battery into the battery compartment
Insufficient massage power	Low battery level	Charge the massage band using the USB Type-C charging cable
	Insufficient battery contact with the device	Reinstall the battery
Incorrect display operation	Short circuit	Contact a certified service center

TECHNICAL SPECIFICATIONS

Name	Head massager
Model	Yamaguchi Headband
Power	3.5 W
Power supply	5V / 700mAh Li-ion battery
Power supply for charging	DC 5 V / 2 A
Type of massage	Vibrating
Vibration frequency	Up to 12000 rpm
Massage areas	Entire head, neck and legs
Number of auto programs	5
Number of massage fingers	11 pcs
Display	Yes
Voice assistant	Joy Assistant
Body material	ABS-plastic
Massage elements material	Silicone
Dimensions (L×W×H)	16×15.2×3 cm
Net weight	120 g
Timer	15 minutes
Noise level	40 dB
Charging time	1.5 hours
Continuous operation time	8.5 hours

YAMAGUCHI

Japan

YAMAGUCHI®

Product of "Yamaguchi Electronics Solutions"
Misumi, Nagato, Yamaguchi, Japan

EAC