SOOCAL | D3 Pro

User Manual

Contents

EN	1
DE	9
FR	18
ES	27
IT	36
RU	45
PL	54
KR	62
JP	70
ZH	78

Thank you for choosing Soocas oral care series: Soocas D3 Pro Electric Toothbrush. Please read this manual carefully and the correct method of usage for Soocas Aura Electric Toothbrush to fully understand its performance to better clean and care for your oral health.

Product Name: Soocas D3 Pro Electric Toothbrush

Model:D3 Pro

Input Power: 5V==0.5A

Charging Time: 6 Hours

Charging Method: Type-C

Product Name: Soocas D3 Pro UVC Toothbrush Sanitizer

Applicable model: D3 Pro

Input Power: 5V==0.5A

Charging Time: 6 hours

Charging Method: Type-C charging

Place of Origin: Shenzhen, China

Manufacturer: Shenzhen Soocas Technology Co., Ltd

Address: Room 1101, 1102, 1103, 1104, Building 2 Chongwen Area,

Nanshan iPark, No. 3370 Liuxian Avenue, Fuguang Community,

Taoyuan Sub-District, Nanshan District, Shenzhen, China

Made in China

All products bearing this symbol are waste electrical and electronic equipment (WEEE as in directive 2012/19/EU) which should not be mixed with unsorted household waste. Instead, you should protect human health and the environment by handing over your waste equipment to a designated collection point for the recycling of waste electrical and electronic equipment, appointed by the government or local authorities.

Correct disposal and recycling will help prevent potential negative consequences to the environment and human health. Please contact installer or local authorities for more information about the location as the well as terms and conditions of such collection points.

User Manual Version: V1.0

Contents

Safety Tips	1
Product Introduction	2
Usage Guide	3
Environmental Protection	7
Troubleshooting	9

Safety Tips

When using this product for the first time, it is common to experience minor discomfort with some gum bleeding. This is because the gums are not quite used to the high-vibration cleaning method. It is recommended to use gentle mode for beginners and expects this discomfort to continue for 1-2 weeks upon initial usage. If bleeding persists for longer than 2 weeks, please stop using the product immediately and consult professional assistance from your dentist.

Warning

- For users with teeth correction, dental repair, dental implants, etc., please use this product under the guidance of a dentist or doctor.
- This product meets the safety standards of electromagnetic devices. For any issues related to pacemakers or other implantable devices, please consult a doctor or the manufacturer of the implant before use.
- This product is only designed for cleaning teeth, gums, and the tongue. Do not use it for any other purpose. If you feel any discomfort or pain, please stop using this product and seek professional assistance from a dentist or a doctor.
- Do not immerse the charging port in water or any liquid.
- Do not plug in the power supply if there is liquid on the charging port or power plug.
- This product has a built-in rechargeable battery. Do not throw it into the fire or charge it in a high-temperature environment.
- Do not knock or drop this product. If the product is damaged due to collision or falling, please stop using it.
- Do not modify, disassemble or repair this product by yourself.
- This device is not suitable for people with physical, sensory, or mental disabilities, and lack of experience and knowledge (including children) to use unless they are supervised or instructed by a person responsible for their safety.
- Children should not play with this product.
- Cleaning and maintenance should not be done by children without supervision.
- Use a safe and compliant power adapter for charging.

Cautious

- Please check for obvious damage before using the product. Do not use it if it is damaged or the machine has suffered a severe fall.
- Make sure that the charging port and charging cable are dry before charging; otherwise, it may cause an electric shock or short circuit.
- This product has a built-in battery, do not throw it into the fire.
- Do not place it in a high-temperature or damp place, and avoid direct sunlight.
- The charging base cannot be used in outdoor direct sunlight.
- Do not use the toothbrush when soaking in a bath or shower.
- Do not share the brush head with others.
- Do not place metal objects on the charging base while it is powered on.
- Note: The product and accessory illustrations in the manual are for reference only. Due to product updates and upgrades, the actual product may differ slightly from the illustration. Please refer to the actual product.
- This product is intended for home-use only. Do not use it in a commercial setting.
- Only use and charge the appliance at temperature between 5°C and 35°C.

Product Introduction

Brush Head
Protecting Cover

D3 Pro
Brush Head

User Manual

Transmission Shaft

Toothbrush Handle

On/Off Button
(Turn On/Off &
Switch Modes)

Gentle Mode

Clean Mode

Brightening Mode

Polishing Mode

Power Indicator

Charging
Cable


Wall Mount
(Placed under
the Brush
Head Box)

USB-C
Charging
Port

Dust-Proof
Rubber Plug

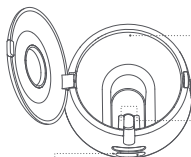Air-Drying
Sterilization
Box

Air Drying
Atmosphere
Lamp (White)
Sterilization
Atmosphere
Lamp (Purple)


Power
Indicator

Dust-Proof
Silicone Plug

Usage Guide

1. Charging the electric toothbrush


Please charge with the matching charging cable before using it for the first time. Built-in rechargeable battery: it takes about 6 hours to fully charge and the specific charging time may differ due to the remaining power and charging environment. The battery temperature must be 0 - 45 °C before charging.

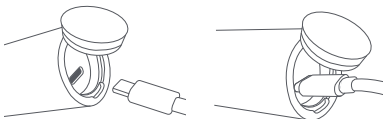Note: In order to ensure good waterproofness and safe use, when the product is not charged, please ensure that the dust-proof rubber plug of the toothbrush has been completely covered at the tail of the handle.

2. Installing the brush head

To begin, remove the brush head and securely attach it. Ensure the bristles are facing forward, then insert the transmission shaft into the hole at the base of the brush head. You'll hear a click when it's properly in place.

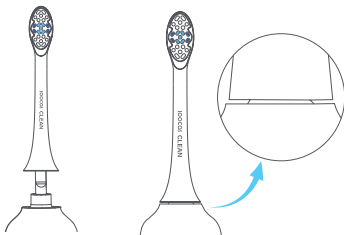3. Modes and functions of toothbrush

(1) Briefly press the on/off button to power up the toothbrush and start brushing your teeth. Within 15 seconds of activation, you can switch between modes by briefly pressing the on/off button again.

(2) In the startup state: Press and hold for 2S to pause. If no action occurs within 30 seconds, the device will shut down, and if the Power Switch is pressed once again within 30 seconds, the device will once again start up and continue the countdown of the last action after the device is switched on.

- | | | |
|---|-------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| ○ | | Gentle Mode: Gentle brushing, caring for new and sensitive users |
| ○ | | Clean Mode: Comfortable and powerful, protects and improves gum health |
| ⌚ | | Brightening Mode: Intensive cleaning power, restores the natural brightness |
| ✦ | | Polishing Mode: High oscillation of the brush head for deep scrubbing of stubborn stains |
| □ | | Battery Reminder Lamp
Red light on: Low power, please charge in duetime
Red or green light flashing: Device is charging
Green light on: The equipment is fully charged |

Zone Change Reminder

Every 30 seconds, it vibrates to pause and to remind replacing the brushing position.

Mode memory function

Every time you turn on it, it is in the same mode as the last time it was used.

Long endurance


The longest battery life of D3 Pro after being fully charged is about 180 days, and the battery life will vary depending on the modes and powers used.

4. Installation and use of sterilization box


(1) Installing the Wall Mount

Attach the wall mount firmly to a suitable surface, including smooth materials like ceramic tiles and glass. Clean and dry the chosen spot on the surface, peel off the backing sticker from the wall mount, align it using the horizontal reference line, and firmly press it onto the surface to ensure a stable attachment. Allow it to set for 24 hours before installing the air-drying sterilization box.

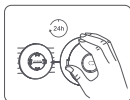
(2) Installing the Air-Drying Sterilization Box

Align the air-drying sterilization box with the wall mount and gently push it into place. Adjust the box's angle as needed.

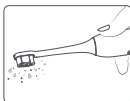
(3) Using the Air-Drying Sterilization Box

After brushing your teeth, gently shake off the residual water droplets on the brush head. Then open the lid of the air-drying sterilization box, hang the bristles inward in the box, close the lid, and the air-drying sterilization box will automatically start air drying for 35 minutes and ultraviolet sterilization for 5 minutes.

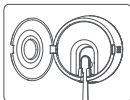(4) Forcibly stop the sterilization box

If you wish to deactivate these functions while in operation, open and close the cover within 5 seconds of starting.

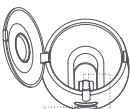Power Indicator Atmosphere Lamp

(5) Indicator light


Working status: The white atmosphere lamp will be on for 35 minutes and then the purple atmosphere lamp will be on for 5 minutes;
Low power: The red power indicator flashes until it is turned off;

Power exhausted: The box cover or opening cover is not started and the red power indicator flashes;
Damage of air drying function: After closing the cover, the white atmosphere lamp flashes quickly for 5 seconds and then goes out;
Damage of sterilization function: After closing the cover, the purple atmosphere lamp flashes quickly for 5 seconds and then goes out.

Note:

- (1) Please note that the air-drying sterilization box is not suitable for washing.
- (2) The wall rack cannot be torn off and pasted for the second time to avoid affecting the firmness after pasting.
- (3) The air-dried sterilization box cannot be paired with other toothbrushes, so as to avoid excessive load-bearing which affects the product life.


5、 Proper Brushing Techniques

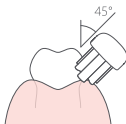- (1) Guide the brush head slowly from tooth to tooth with a 40-degree angle between the bristles and the gum lines.

- (2) Press the On/Off button for the preferred brushing mode.

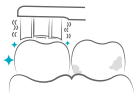- (3) The toothbrush will vibrate to indicate the user has to brush a new area of the mouth.

- (4) Keep the brush head to the teeth closely without bending the bristles.

Environmental Protection

When the product is at the end of its service life, according to local laws and regulations, one must remove the internal battery and undergo proper recycling and disposal procedures before the parts can be disposed of. If you are having trouble displacing the battery, please dispose of the product in accordance with the local waste clarification regulation.

Warning

- Do not throw batteries into the fire to prevent battery explosion or the release of toxic and harmful substances.
- Do not short-circuit the battery to prevent battery combustion and fire.
- Unless the intention is to discard the product, do not disassemble the product parts to prevent damaging the product.
- Strict safety procedures must be taken place. Before removing the battery, the power connection must be disconnected.
- Before discarding the product, the battery must be removed from the product.
- The product must be powered off when removing the battery.
- Batteries should be disposed of safely.

Remove the rechargeable battery

1. Use a flathead screwdriver to pry open the handle end cap.
2. Detach the cover and use a flathead screwdriver to pray open the inner motor. Pop out the motor.
3. Cut the red wires on the inner motor.


Remove the rechargeable battery:


1. Insert a straight screwdriver into the gap between the tail cover and the main body of the handle, and pry up until the left and right buckle positions are opened at the tail cover.


2. Use one end of the shaft against the table and apply force downward.

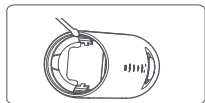
3. While conducting the second step, pry out the two symmetrical buckles of the inner bracket with a straight screwdriver, and push out the movement with downward force on the body.

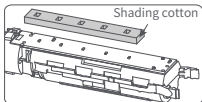4. Remove the shading cotton.

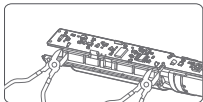
5. Forcibly cut between the plastic and PCB indicated by the arrow with bevel pliers until the solder pin of the battery is cut.

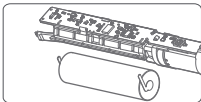
6. Take out the battery and put it in the designated recycling place.


How to remove the battery of the sterilization box:

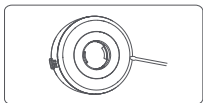1. Insert a straight screwdriver into the gap between the decorating part of the bottom shell and the bottom shell, and pry up until the decorating part of the bottom shell is opened and removed.

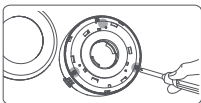2. Remove the decorative ring from the bottom shell and remove the three screws with a Phillips screwdriver.

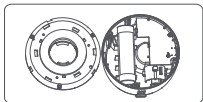3. After removing the three screws, open the back shell.

4. Cut the two power cords with cutting pliers, remove the battery and place it in the specified recycling place.

Troubleshooting

Problem	Possible Causes	Solutions
Unable to boot	Disabled while charging	Stop charging
Unable to boot, red light flashes 3 times	Low power	Charge in time
Unable to boot, red light flashes 6 times	The temperature of the battery is too high or too low	Please use at suitable temperature (-20°C - 60 °C)
Failing to charge	Abnormal temperature of battery	Please charge at appropriate temperature (0 - 45°C)
Unable to install the brush head	There is a foreign matter in the hole of the brush head	Remove the foreign matter in the hole of the brush head
	There is a viscous foreign matter on the transmission shaft	Remove viscous foreign matter
	It is not the brush head in the correct version	Please use the brush head in the correct version
	The insertion direction is reversed	Please refer to the brush head installation instructions
Weak vibration of the brush head	The brush head is not installed correctly	Reinstall the brush head
Abnormal/unlit indicator light	Battery depletion	Charge in time
	Damage in circuit board	Return to after-sales for maintenance
Abnormal/unlit indicator light of the air-dried sterilization box	Insufficient power	Charge in time
	The white atmosphere light flashes quickly and then goes out	Return to after-sales for maintenance
	The purple atmosphere light flashes quickly and then goes out	Return to after-sales for maintenance

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