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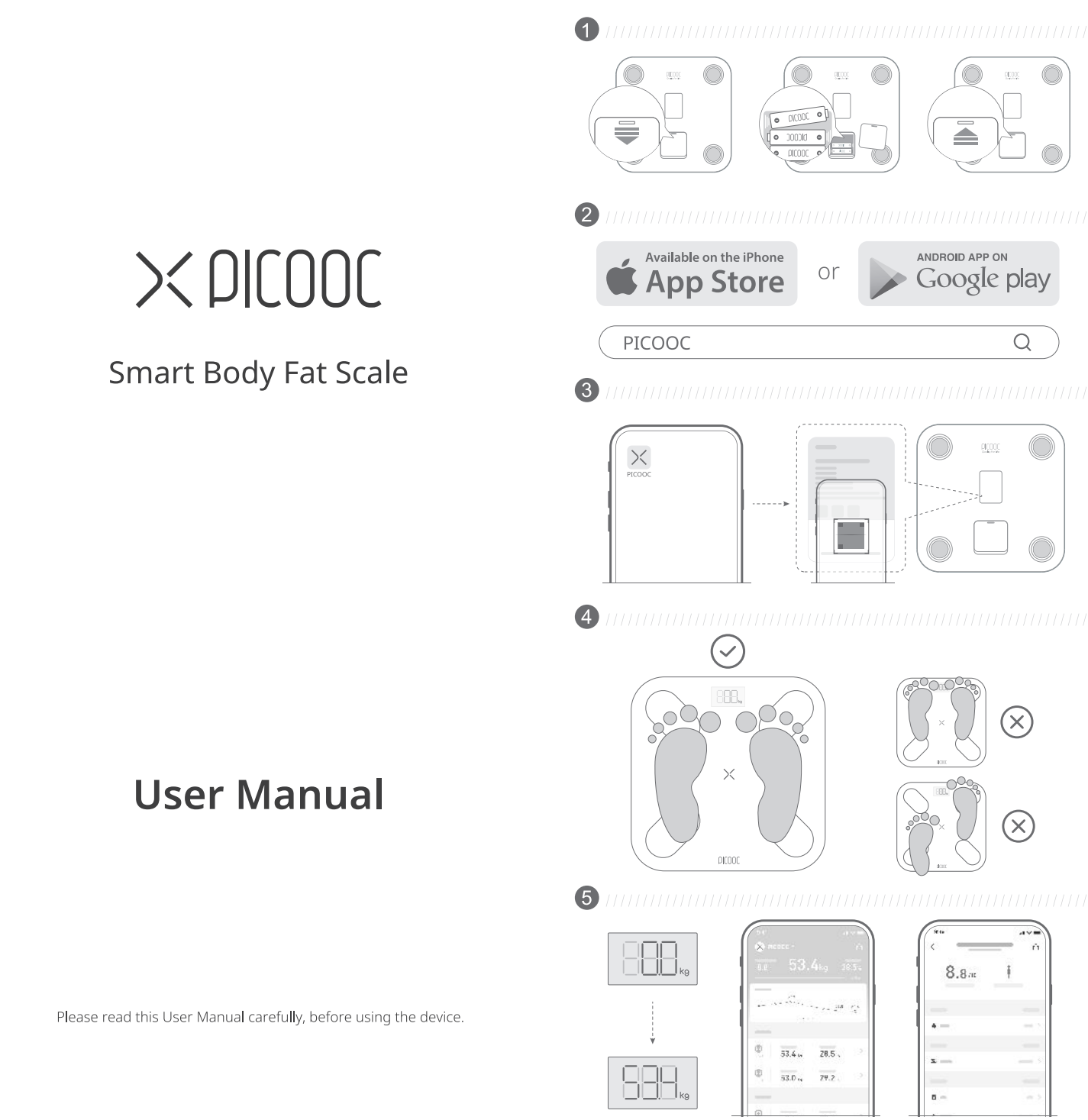
English

Security Information

- ⚠ Patients with pacemakers or other medical equipment should not use this device!
- ⚠ If you are pregnant, please consult a doctor for advice before using!
- ⚠ Do not place the device on slippery surfaces when using!
- ⚠ The device is covered by glass. In order to prevent injury, please make sure that the glass is not damaged before using.
- ⚠ Any data or advices in the APP are for reference only, and they do not replace or stand for professional diagnosis of medical institutions.
- ⚠ Please keep the manual properly, and offer it to others before they use the device, the manual should be included for sales.

For more security considerations, Please check your local dealer supporting websites or visit <https://www.picooc.com>

Measuring Instructions

Measurement error may be larger for the following users:

- Growing children
- Elderly people
- Fever patients due to cold and other reasons
- Osteoporotic patients with very low bone mineral density
- Edema patients
- Artificial dialysis patients
- People who are involved with fitness or sports (or similar activities)

Please avoid using the scale in the following situations, as far as possible:

- After strenuous exercise
- After sauna or a bath
- After heavy drinking
- After ingesting significant food and/or water

Other measurements

Light object weighing

Pressing the scale surface until the LED screen display 0.00, place the object on the scale and the LED screen would display its weight (measurement range: 200 g to 10 kg)

Heart rate measurement

Keep measuring posture for a few seconds after measuring body fat, LED screen will show resting heart rate value (measurement range: 40-125 beats/min)

Balance ability test

Please follow the operation guide and particularly be sure to comply with the safety precautions in the app

Warranty and Support

You may consult your dealer for warranty and service support policies of the device.

You may check the dealer information on supporting websites and learn about more information of the device. (<http://www.picooc.com>)

Deutsch

Please read this User Manual carefully, before using the device.

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