
75mm

210mm

DR.BEI Portable Water Flosser GF1
User Manual
Please read this instruction manual carefully before using the product and save it for future use.

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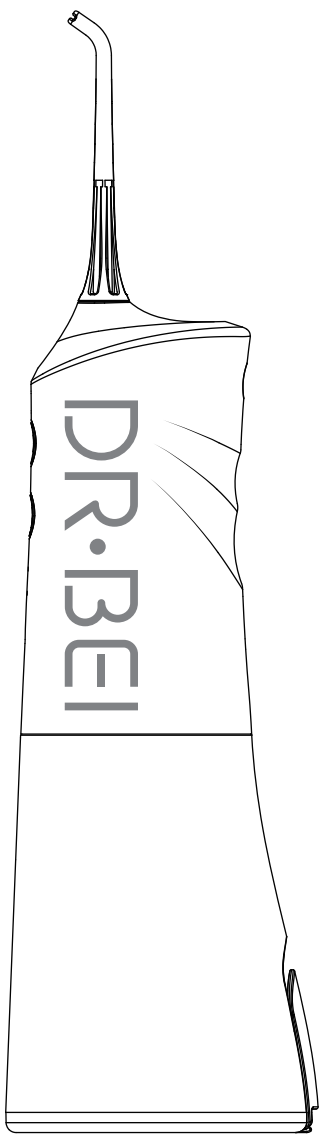IMPORTANT SAFE GUARDS
All information in this section is provided to instruct you in the correct and safe operation of this device so that you can avoid injuries to yourself and others as well as damage to property.


READ ALL INSTRUCTIONS BEFORE USING

DANGERS

- Toddlers, infants, people who are unable to operate should not allowed use this appliance. Please keep it out of reach of children.
- Never alter the appliance. Also do not try to disassemble or repair it.
- Do not place or store the product where it can fall or be pulled into tub or sink.
- Do not reach for the product that has fallen into water. Unplug immediately.
- Do not use the product while bathing.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- This appliance contains batteries that are non-replaceable.

WARNINGS

- Oral irrigators can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.
- Children shall not play with the appliance.
- The appliance is only to be used with the power supply unprovided with the appliance.
- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery.
- The battery is to be disposed of safely.
- Do not fill the tank with water warmer over 40°C.
- Never use the adaptor if it has a damaged cord or plug.
- Keep the cord away from heated surfaces.
- Please do not use it if you lose the nozzle.
- Do not submerge the device in water.
- Do not use for any purpose other than oral cleaning.
- Keep nozzle away from small children. The nozzle is small enough to be swallowed and can be become lodged in the throat.
- Certain types of mouthwash may damage the equipment by causing the housing or water reservoir to crack. Therefore, do not use mouthwash in the water reservoir.
- People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
- Do not press the nozzle too firmly against your teeth or gums.
- Doing so may injure your gums.

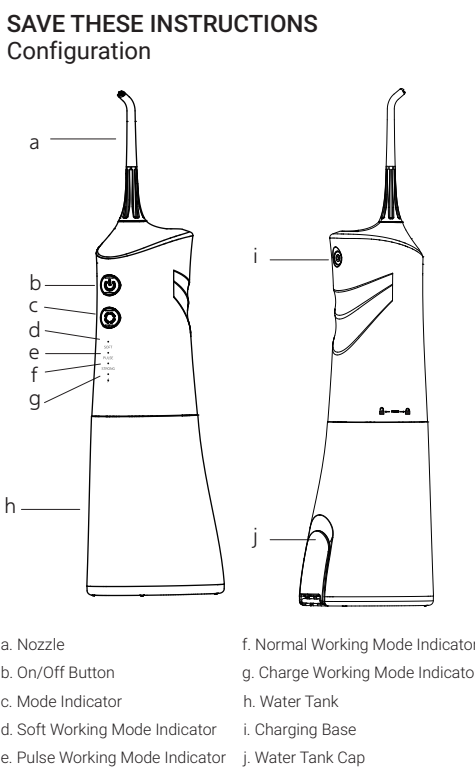
Note:
There may be some liquid remaining inside the main appliance before use (tank, water suction tube, etc), this is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.

How to use
Charging the device (Fig. 1)

1. Take the charger base cover down and insert the appliance plug securely into the appliance socket (make sure no liquid remaining inside the socket before insert the appliance plug).
2. Insert the charging adaptor securely into the power outlet.
3. The constant flash a red LED indicates charging is in process.
4. The constant red LED will die out in 2 minutes indicates that the battery is fully charged.

Note:
During charging in progress, the indicator is lit in red. After fully charged the indicator is lit in blue. If the indicator is blinking red during use, it means battery exhausted, please recharge the appliance for continuous use. It takes 5.5 hours until fully charged. Make sure there is no liquid remaining in the electrical outlet before plugging.

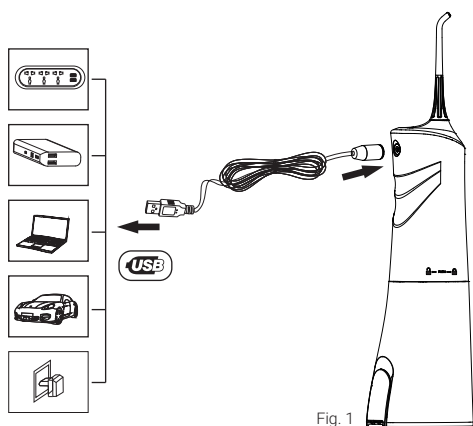Use the water flosser

1. Open the water tank cap to fill the tank with water which temperature not exceed 40°C, then close the cap. (Fig. 2)
2. When first use please add small amount of water to water outlet. (Fig. 3)
3. Attach nozzle onto the unit, rotate 45°. (Fig. 4)
4. How to put on the tank. (Fig. 5)

After using the irrigator (Fig. 10)

1. Open the water tank cap and dispose of the water in the tank.
2. Press "On/Off" key to turn on the unit, and drain the water off completely.
3. Press "On/Off" key again to turn off the unit, and close the water tank cap.
4. Wipe the device with dry cloth.

Specifications

Product Name	DR.BEI Portable Water Flosser GF1
Model	DR.BEI GF1
Rated Voltage	3.7V rrr
Rated Power	5W
IP Rating	IPX7
Charging Time	Full recharge takes at least 5 hours
Rated Input	5V rrr 1A
Place of Origin	China

6. First press "Mode" key to select the mode you desired, then the on key. During irrigating, you can change the working mode (Fig. 7)

SOFT: For gentle cleaning of sensitive gums. The second indicator is white.

PULSE: For cleaning and massaging the gums. The third indicator is white.

STRONG: For normal cleaning. The first indicator is white.


7. Hold the unit in an upright position with nozzle tip point toward your teeth.

8. Keep your mouth opened slightly so that the water can flow out.

9. Press "On/Off" key to start irrigating. (Fig. 8)

10. Attach nozzle onto the unit, rotate 45°. (Fig. 9)

11. Once irrigating is completed, press "On/Off" key again to turn off the device, press the button.

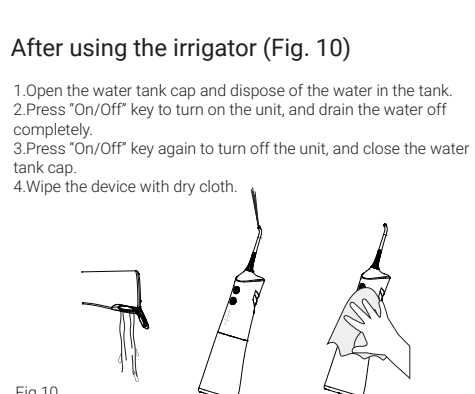After-sales Service
If the product has any fault mentioned in the List of Product Faults, please contact a local dealer or an after-sales service provider.

Warranty restrictions (the warranty does not cover the following):

1. Unauthorized repair, misuse, damage by shock, neglect, abuse, permeation with liquid, accident, or alteration, or incorrect use of any part not for the product, or torn or altered labels or anti-counterfeit marks.
2. Expiry of warranty.
3. Damage due to force majeure.
4. Inconsistency with the faults mentioned in the List of Product Faults.
5. Faults of the product or its part which is resulted from deliberate damage and not mentioned in the List of Product Faults.

Note:
To prevent bacteria propagating, make sure there is no remaining water inside the unit. If you're not going to use it for a long time, be sure wipe and dry the appliance ASAP before store it.

Maintenance
Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents. (e.g. vinegar or lime scale removers), since they may damage the device. Do not wash using hot water of over 40°C. Keep the device away from high temperature area and avoid direct sunlight.

Main unit

1. Remove and the dirt clean with a liquid soap and wipe off using a dry cloth.
2. Do not submerge the main unit in water to clean.
3. Dry after washing with water.

Water tank cap
Wipe off any grime and dirt with a cloth.

Nozzle

1. Wash with water and wipe with soft cloth.
2. Replace the nozzle every 6 months.

Water suction tube

1. Wash with water and wipe with soft cloth.
2. Do not bend, pull or twist the water suction tube.

Tank

1. Wash with water and wipe with soft cloth.
2. Dry the tank when not using for a long period.

Troubleshooting

Problem

1. It does not work.
2. The battery has reached the end of its life cycle.
3. Water pressure is very low.
4. Water does not come out.

Possible Causes

1. You've just purchased the product or it has been left unused for more than 3 months.
2. The battery has reached the end of its life cycle.
3. There is dirt stuck in the nozzle.
4. The nozzle is deformed.
5. There is dirt stuck in the filter.
7. The tank is empty.
8. The irrigator is lifted too much while in use.
9. Valve plate slight deformation.

Solution

1. Charge the unit for at least 5 hours.
2. Contact service center for support.
3. Replace the nozzle.
4. Clean the filter.
5. Fill up tank.
6. Use with the appliance in a vertical state.
7. Add small amount of water to water outlet.

List of Product Faults

Name	Faults
DR.BEI Portable Water Flosser	The product fails to switch on or off. No water comes out when the water tank has water and there is no dirt stuck in the nozzle. The power indicator is not functioning when the battery is being charged. The battery does not charge properly.

After-sales Service: service@doctorbei.com
Manufacturer: Wuxi Qinghe Xiebei Technology Co., Ltd.
Address: Room 1106, 11th Floor, NO.2 of Building 530, Binhu District, Wuxi City, Jiangsu, P.R. China

Environmental Protection
When the product has reached the end of its service life, it may be necessary to remove the internal battery for recycling according to relevant national or local laws and regulations. If you cannot take out the battery, please hand over the product to the recycling center designated by the government.

CE EU Declaration of Conformity
We Wuxi Qinghe Xiebei Technology Co., Ltd., hereby, declares that this equipment is in compliance with the applicable Directives and European Norms, and amendments.

WEEE Disposal and Recycling
Correct Disposal of this product. This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

FC Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/ TV technician for help.

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